



**Glidewell HT™**  
IMPLANT SYSTEM

Product Catalog

[glidewell-ht.com](http://glidewell-ht.com)

|           |                                               |           |                                            |
|-----------|-----------------------------------------------|-----------|--------------------------------------------|
| 5         | Technical Specifications                      | 22        | Drivers & Accessories                      |
| 6         | Implants                                      | <b>23</b> | <b>Guided Surgical Instrumentation Kit</b> |
| 7         | Healing Components                            | 24        | Guided Surgical Kit Overview               |
| 8         | Impression Components                         | 26        | Guided Drills                              |
| 10        | Stock Abutments                               | 28        | Guided Taps & Accessories                  |
| 13        | Multi-Unit Abutments                          | <b>30</b> | <b>Prosthetic Instrumentation Kit</b>      |
| 15        | Multi-Unit Accessories                        | 31        | Prosthetic Instrumentation Kit Overview    |
| 16        | LOCATOR® Abutment System                      | <b>32</b> | <b>Product Packaging</b>                   |
| <b>18</b> | <b>Standard Surgical Instrumentation Kits</b> | <b>33</b> | <b>Policies and Warranty</b>               |
| 19        | Surgical Kit Overview                         | <b>34</b> | <b>Parts Index</b>                         |
| 20        | Drills & Taps                                 |           |                                            |

***Glidewell HT™ Implant System is a trademark of Prismahtik Dentalcraft, Inc.***

Copyright © 2025, Prismahtik Dentalcraft, Inc. Prismahtik Dentalcraft, Inc. is not responsible for any damages or other liabilities (including attorney fees) resulting, or claimed to result in whole or in part, from actual or alleged problems arising out of the use of this information. The techniques, procedures and theories presented herein are provided in good faith and believed to be correct as of the date hereof. Any dental professional viewing this presentation must make his or her own decisions about the use of the materials and techniques for specific situations.

*No representations as to the completeness or accuracy of this information is given, and no representations or warranties, either expressed or implied, of merchantability, fitness for a particular purpose or of any other nature are made here under with respect to the information or the product to which information refers.*

*All third-party trademarks are property of their respective owners.*

# The Glidewell HT™ Implant System

The **Glidewell HT™ Implant System**, formerly known as the Hahn™ Tapered Implant System, is a premium implant solution that simplifies surgery and provides you with unrivaled support from the most experienced dental lab in the U.S. With a 99.2% long-term success rate backed by an implant-to-crown lifetime warranty, the Glidewell HT Implant promotes success in implant dentistry while lowering your surgical and restorative costs.<sup>1</sup>

- **Simple and efficient** — Easy to use, with a streamlined surgical protocol, length-specific drills, and new-and-improved surgical kits based on eight years of clinical feedback
- **Cut your costs** — Priced at a fraction of comparable implant systems — and saves you 20% on your lab bill when you restore your implant case with Glidewell\*
- **Clinically proven** — 99.2% success rate and 0.2 mm mean bone loss<sup>1</sup>
- **High primary stability** — Deep, sharp threads maximize initial stability and engage bone where directed

“The **Glidewell HT™ Implant System** is engineered to help dentists provide implant treatment for more patients through **ease of use, reduced costs** and our **unwavering commitment to support your practice** — from implant placement to final restoration.”

**Jim Glidewell, CDT**  
Founder and President of Glidewell



## About the Manufacturer

Prismatik Dentalcraft was established in 2006 and includes a carefully assembled team of experts with a proven track record in the design, engineering and manufacture of dental implants. Bolstered by a support staff of highly respected researchers, material scientists, clinical specialists, and dental technicians, Prismatik is dedicated to advancing implant therapies by combining proven treatment protocols with progressive materials, technologies and techniques.



### Vertical Integration

Our ownership of the entire manufacturing process behind our implant products ensures quality and helps reduce costs for our customers.



### State-of-the-Art Equipment

Our Swiss-type lathes and CNC milling machines are ideal for implants and prosthetics requiring extreme precision.



### Made in the USA

Our ISO-certified facility in Irvine, California, operates under FDA Current Good Manufacturing Practices (CGMPs).

\*Discount offered only at Glidewell and cannot be combined with any other special offers. Case must include an implant-level or multi-unit abutment-level impression with a Glidewell HT transfer coping or a digital scan with a Glidewell HT scan body. Impressions over cementable abutments are not eligible for discount.

1. Kerr M, Allen B, Park N. Clinical and radiographic evaluation of tapered implants with an aggressive reverse buttress thread and crestal microthreads: a retrospective study. For the full report, visit [glidewell.com/ht-2-year](http://glidewell.com/ht-2-year).

## PROVEN, TIME-TESTED IMPLANT DESIGN FEATURES

# 99.2%

SUCCESS RATE<sup>1</sup>

# 0.2 mm

MEAN BONE LOSS IN 2-YEAR STUDY<sup>1</sup>

### REDUCED MARGINAL BONE STRESS<sup>2</sup>

Convergent implant neck design minimizes forces on cortical bone.

### ENHANCED BIOLOGICAL SEAL<sup>2,13</sup>

Platform-shifted abutment interface and conical connection provide excellent seal, stability and tissue preservation.

### RBM SURFACE TECHNOLOGY<sup>3-5</sup>

Resorbable blast media (RBM) surface treatment results in textured, ultra-clean surface that promotes osseointegration.

### SWIFT INSERTION<sup>6</sup>

Dual-lead thread pattern with self-tapping grooves.



### CRESTAL BONE AND SOFT-TISSUE PRESERVATION<sup>2,7-9</sup>

Machined collar surface leads to 25-times less biofilm attachment and mitigates peri-implantitis. Coronal microthreads reduce crestal bone loss.

### HIGH PRIMARY STABILITY AND DIRECTIONAL CONTROL DURING PLACEMENT<sup>10,11</sup>

Deep, sharp threads maximize initial stability for immediate load opportunities and long-term success.

### TAPERED IMPLANT BODY<sup>12</sup>

Tapered body shape — a design feature originally conceived by Dr. Jack Hahn — eases placement in anatomically constricted areas.

# Technical Specifications



Ø3.0 mm

11.5 mm  
13 mm  
16 mm

Ø3.5 mm

8 mm  
10 mm  
11.5 mm  
13 mm  
16 mm

Ø4.3 mm

8 mm  
10 mm  
11.5 mm  
13 mm  
16 mm

Ø5.0 mm


8 mm  
10 mm  
11.5 mm  
13 mm  
16 mm

Ø6.0 mm


8 mm  
10 mm  
11.5 mm

Ø7.0 mm

8 mm  
10 mm  
11.5 mm

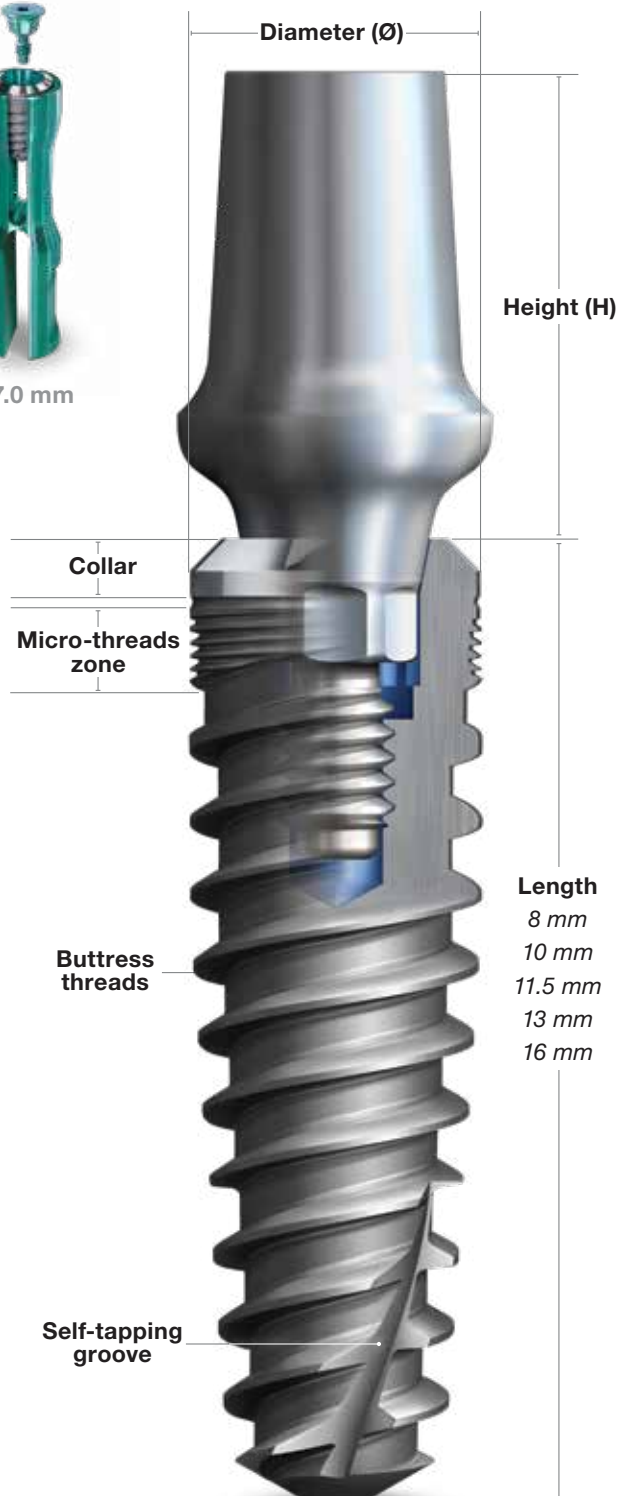

Glidewell HT™ Implants allow for precise control during placement, engage a maximum amount of bone, and achieve a high degree of primary stability in a wide variety of clinical situations. The unique thread pattern with self-tapping grooves is sharply angled for swift insertion. All Glidewell HT Implants feature a conical internal hex connection for a secure prosthetic seal. Matching components are color-coded by implant diameter and/or restorative platform. To help minimize the overall number of system components, 3.5 mm and 4.3 mm diameter implants share the same restorative platform.

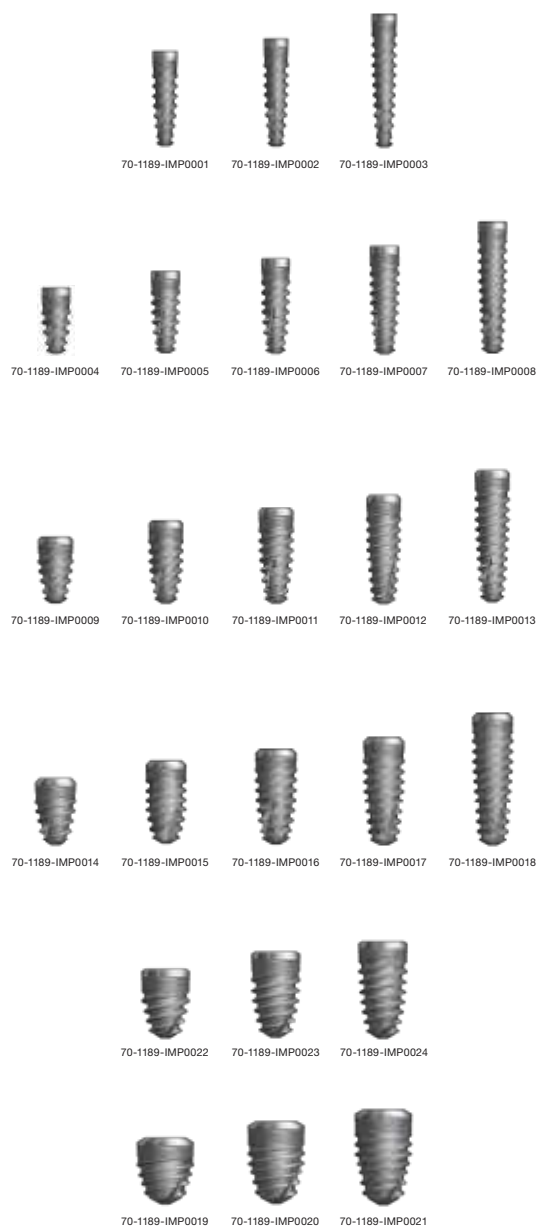


| Description                  | Code            |
|------------------------------|-----------------|
| <b>Glidewell HT™ Implant</b> |                 |
| Ø3.0 x 11.5 mm               | 70-1189-IMP0001 |
| Ø3.0 x 13 mm                 | 70-1189-IMP0002 |
| Ø3.0 x 16 mm                 | 70-1189-IMP0003 |
| Ø3.5 x 8 mm                  | 70-1189-IMP0004 |
| Ø3.5 x 10 mm                 | 70-1189-IMP0005 |
| Ø3.5 x 11.5 mm               | 70-1189-IMP0006 |
| Ø3.5 x 13 mm                 | 70-1189-IMP0007 |
| Ø3.5 x 16 mm                 | 70-1189-IMP0008 |
| Ø4.3 x 8 mm                  | 70-1189-IMP0009 |
| Ø4.3 x 10 mm                 | 70-1189-IMP0010 |
| Ø4.3 x 11.5 mm               | 70-1189-IMP0011 |
| Ø4.3 x 13 mm                 | 70-1189-IMP0012 |
| Ø4.3 x 16 mm                 | 70-1189-IMP0013 |
| Ø5.0 x 8 mm                  | 70-1189-IMP0014 |
| Ø5.0 x 10 mm                 | 70-1189-IMP0015 |
| Ø5.0 x 11.5 mm               | 70-1189-IMP0016 |
| Ø5.0 x 13 mm                 | 70-1189-IMP0017 |
| Ø5.0 x 16 mm                 | 70-1189-IMP0018 |
| Ø6.0 x 8 mm                  | 70-1189-IMP0022 |
| Ø6.0 x 10 mm                 | 70-1189-IMP0023 |
| Ø6.0 x 11.5 mm               | 70-1189-IMP0024 |
| Ø7.0 x 8 mm                  | 70-1189-IMP0019 |
| Ø7.0 x 10 mm                 | 70-1189-IMP0020 |
| Ø7.0 x 11.5 mm               | 70-1189-IMP0021 |

(Implants include Titanium Cover Screw.)

## Description

## Code

### Glidewell HT™ Implant Healing Abutment

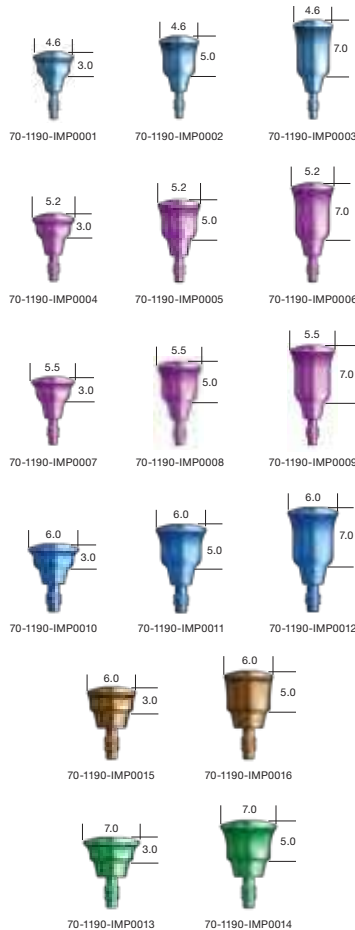

|                      |                 |
|----------------------|-----------------|
| 3 mmH - Ø3.0 Implant | 70-1190-IMP0001 |
| 5 mmH - Ø3.0 Implant | 70-1190-IMP0002 |
| 7 mmH - Ø3.0 Implant | 70-1190-IMP0003 |
| 3 mmH - Ø3.5 Implant | 70-1190-IMP0004 |
| 5 mmH - Ø3.5 Implant | 70-1190-IMP0005 |
| 7 mmH - Ø3.5 Implant | 70-1190-IMP0006 |
| 3 mmH - Ø4.3 Implant | 70-1190-IMP0007 |
| 5 mmH - Ø4.3 Implant | 70-1190-IMP0008 |
| 7 mmH - Ø4.3 Implant | 70-1190-IMP0009 |
| 3 mmH - Ø5.0 Implant | 70-1190-IMP0010 |
| 5 mmH - Ø5.0 Implant | 70-1190-IMP0011 |
| 7 mmH - Ø5.0 Implant | 70-1190-IMP0012 |
| 3 mmH - Ø6.0 Implant | 70-1190-IMP0015 |
| 5 mmH - Ø6.0 Implant | 70-1190-IMP0016 |
| 3 mmH - Ø7.0 Implant | 70-1190-IMP0013 |
| 5 mmH - Ø7.0 Implant | 70-1190-IMP0014 |

(Healing Abutments require a Glidewell HT™ Implant Prosthetic Driver for placement; please see page 22.)

### Glidewell HT™ Implant Temporary Abutment



|                  |                 |
|------------------|-----------------|
| Ø3.0 Implant     | 70-1190-PRA0025 |
| Ø3.5/4.3 Implant | 70-1190-PRA0026 |
| Ø5.0 Implant     | 70-1190-PRA0027 |
| Ø6.0 Implant     | 70-1190-PRA0062 |
| Ø7.0 Implant     | 70-1190-PRA0028 |

(Temporary Abutments include Titanium Screw.)

## Description

## Code


### Glidewell HT™ Implant Non-Engaging Temporary Abutment



|                  |                 |
|------------------|-----------------|
| Ø3.0 Implant     | 70-1190-PRA0021 |
| Ø3.5/4.3 Implant | 70-1190-PRA0022 |
| Ø5.0 Implant     | 70-1190-PRA0023 |
| Ø6.0 Implant     | 70-1190-PRA0063 |
| Ø7.0 Implant     | 70-1190-PRA0024 |

(Temporary Abutments include Titanium Screw.)

### Glidewell HT™ Implant Closed-Tray Impression Coping

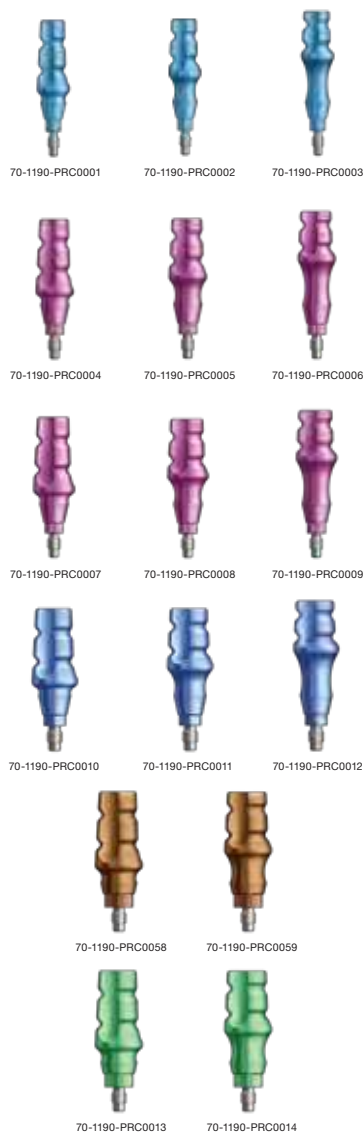

|                      |                 |
|----------------------|-----------------|
| 3 mmH - Ø3.0 Implant | 70-1190-PRC0001 |
| 5 mmH - Ø3.0 Implant | 70-1190-PRC0002 |
| 7 mmH - Ø3.0 Implant | 70-1190-PRC0003 |
| 3 mmH - Ø3.5 Implant | 70-1190-PRC0004 |
| 5 mmH - Ø3.5 Implant | 70-1190-PRC0005 |
| 7 mmH - Ø3.5 Implant | 70-1190-PRC0006 |
| 3 mmH - Ø4.3 Implant | 70-1190-PRC0007 |
| 5 mmH - Ø4.3 Implant | 70-1190-PRC0008 |
| 7 mmH - Ø4.3 Implant | 70-1190-PRC0009 |
| 3 mmH - Ø5.0 Implant | 70-1190-PRC0010 |
| 5 mmH - Ø5.0 Implant | 70-1190-PRC0011 |
| 7 mmH - Ø5.0 Implant | 70-1190-PRC0012 |
| 3 mmH - Ø6.0 Implant | 70-1190-PRC0058 |
| 5 mmH - Ø6.0 Implant | 70-1190-PRC0059 |
| 3 mmH - Ø7.0 Implant | 70-1190-PRC0013 |
| 5 mmH - Ø7.0 Implant | 70-1190-PRC0014 |

(Closed-Tray Impression Copings include Closed-Tray Screw.)


## Description

## Code


### Glidewell HT™ Implant Open-Tray Impression Coping




|                      |                 |
|----------------------|-----------------|
| 3 mmH - Ø3.0 Implant | 70-1190-PRC0015 |
| 5 mmH - Ø3.0 Implant | 70-1190-PRC0016 |
| 7 mmH - Ø3.0 Implant | 70-1190-PRC0017 |

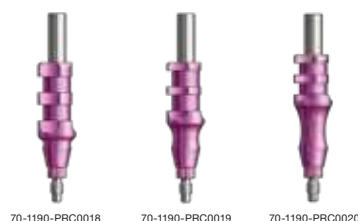


|                      |                 |
|----------------------|-----------------|
| 3 mmH - Ø3.5 Implant | 70-1190-PRC0018 |
| 5 mmH - Ø3.5 Implant | 70-1190-PRC0019 |
| 7 mmH - Ø3.5 Implant | 70-1190-PRC0020 |

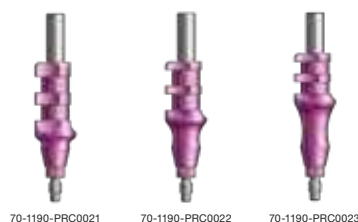

|                      |                 |
|----------------------|-----------------|
| 3 mmH - Ø4.3 Implant | 70-1190-PRC0021 |
| 5 mmH - Ø4.3 Implant | 70-1190-PRC0022 |
| 7 mmH - Ø4.3 Implant | 70-1190-PRC0023 |

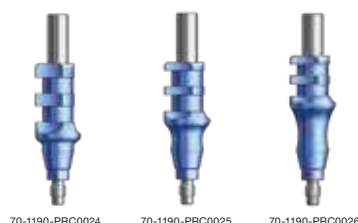

|                      |                 |
|----------------------|-----------------|
| 3 mmH - Ø5.0 Implant | 70-1190-PRC0024 |
| 5 mmH - Ø5.0 Implant | 70-1190-PRC0025 |
| 7 mmH - Ø5.0 Implant | 70-1190-PRC0026 |

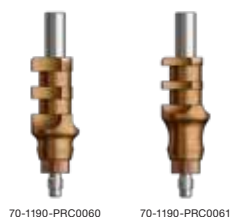


|                      |                 |
|----------------------|-----------------|
| 3 mmH - Ø6.0 Implant | 70-1190-PRC0060 |
| 5 mmH - Ø6.0 Implant | 70-1190-PRC0061 |



|                      |                 |
|----------------------|-----------------|
| 3 mmH - Ø7.0 Implant | 70-1190-PRC0027 |
| 5 mmH - Ø7.0 Implant | 70-1190-PRC0028 |

(Open-Tray Impression Copings include Open-Tray Screw.)

## Impression Components

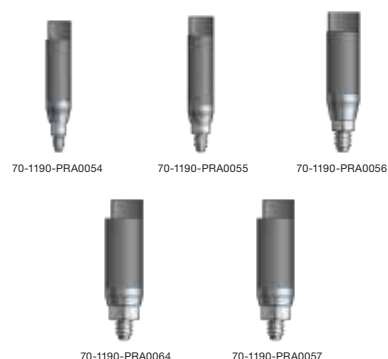

| Description                                     | Code            |
|-------------------------------------------------|-----------------|
| <b>Glidewell HT™ Implant Titanium Scan Body</b> |                 |
| Ø3.0 Implant                                    | 70-1190-PRA0054 |
| Ø3.5/4.3 Implant                                | 70-1190-PRA0055 |
| Ø5.0 Implant                                    | 70-1190-PRA0056 |
| Ø6.0 Implant                                    | 70-1190-PRA0064 |
| Ø7.0 Implant                                    | 70-1190-PRA0057 |

### Glidewell HT™ Implant Analog



|              |                 |
|--------------|-----------------|
| Ø3.0 Implant | 70-1190-PRC0035 |
| Ø3.5 Implant | 70-1190-PRC0036 |
| Ø4.3 Implant | 70-1190-PRC0037 |
| Ø5.0 Implant | 70-1190-PRC0038 |
| Ø6.0 Implant | 70-1190-PRC0062 |
| Ø7.0 Implant | 70-1190-PRC0039 |

## Stock Abutments

### Glidewell HT™ Implant Anterior Esthetic Abutment



|                  |                 |
|------------------|-----------------|
| Ø3.0 Implant     | 70-1190-PRA0037 |
| Ø3.5/4.3 Implant | 70-1190-PRA0038 |
| Ø5.0 Implant     | 70-1190-PRA0039 |

### Glidewell HT™ Implant Posterior Esthetic Abutment



|              |                 |
|--------------|-----------------|
| Ø5.0 Implant | 70-1190-PRA0040 |
| Ø6.0 Implant | 70-1190-PRA0075 |
| Ø7.0 Implant | 70-1190-PRA0041 |

(Esthetic Abutments include Titanium Screw.)

## Description

## Code

### Glidewell HT™ Implant 15° Anterior Esthetic Abutment



|                  |                 |
|------------------|-----------------|
| Ø3.0 Implant     | 70-1190-PRA0042 |
| Ø3.5/4.3 Implant | 70-1190-PRA0043 |
| Ø5.0 Implant     | 70-1190-PRA0044 |

### Glidewell HT™ Implant 15° Posterior Esthetic Abutment



|              |                 |
|--------------|-----------------|
| Ø5.0 Implant | 70-1190-PRA0045 |
|--------------|-----------------|

(Esthetic Abutments include Titanium Screw.)

### Glidewell HT™ Implant Non-Engaging Titanium Abutment



|                            |                 |
|----------------------------|-----------------|
| 4.5 mmH - Ø3.0 Implant     | 70-1190-PRA0001 |
| 4.5 mmH - Ø3.5/4.3 Implant | 70-1190-PRA0002 |
| 4.5 mmH - Ø5.0 Implant     | 70-1190-PRA0003 |



|                        |                 |
|------------------------|-----------------|
| 4.5 mmH - Ø6.0 Implant | 70-1190-PRA0065 |
| 4.5 mmH - Ø7.0 Implant | 70-1190-PRA0004 |

### Glidewell HT™ Implant Titanium Abutment



|                            |                 |
|----------------------------|-----------------|
| 4.5 mmH - Ø3.0 Implant     | 70-1190-PRA0005 |
| 4.5 mmH - Ø3.5/4.3 Implant | 70-1190-PRA0006 |
| 4.5 mmH - Ø5.0 Implant     | 70-1190-PRA0007 |



|                        |                 |
|------------------------|-----------------|
| 4.5 mmH - Ø6.0 Implant | 70-1190-PRA0066 |
| 4.5 mmH - Ø7.0 Implant | 70-1190-PRA0008 |



|                          |                 |
|--------------------------|-----------------|
| 6 mmH - Ø3.0 Implant     | 70-1190-PRA0009 |
| 6 mmH - Ø3.5/4.3 Implant | 70-1190-PRA0010 |
| 6 mmH - Ø5.0 Implant     | 70-1190-PRA0011 |



|                      |                 |
|----------------------|-----------------|
| 6 mmH - Ø6.0 Implant | 70-1190-PRA0067 |
| 6 mmH - Ø7.0 Implant | 70-1190-PRA0012 |

(Titanium Screw sold separately; please see page 13.)

| Description                                                                         |                  | Code            |
|-------------------------------------------------------------------------------------|------------------|-----------------|
| <b>Glidewell HT™ Implant UCLA Gold</b>                                              |                  |                 |
|    | Ø3.0 Implant     | 70-1190-PRA0046 |
|    | Ø3.5/4.3 Implant | 70-1190-PRA0047 |
|    | Ø5.0 Implant     | 70-1190-PRA0048 |
|    | Ø6.0 Implant     | 70-1190-PRA0070 |
|    | Ø7.0 Implant     | 70-1190-PRA0049 |
| <b>Glidewell HT™ Implant Non-Engaging UCLA Gold</b>                                 |                  |                 |
|    | Ø3.0 Implant     | 70-1190-PRA0050 |
|    | Ø3.5/4.3 Implant | 70-1190-PRA0051 |
|    | Ø5.0 Implant     | 70-1190-PRA0052 |
|   | Ø6.0 Implant     | 70-1190-PRA0073 |
|   | Ø7.0 Implant     | 70-1190-PRA0053 |
| <b>Glidewell HT™ Implant UCLA Plastic</b>                                           |                  |                 |
|  | Ø3.0 Implant     | 70-1190-PRA0029 |
|  | Ø3.5/4.3 Implant | 70-1190-PRA0030 |
|  | Ø5.0 Implant     | 70-1190-PRA0031 |
|  | Ø6.0 Implant     | 70-1190-PRA0074 |
|  | Ø7.0 Implant     | 70-1190-PRA0032 |
| <b>Glidewell HT™ Implant Non-Engaging UCLA Plastic</b>                              |                  |                 |
|  | Ø3.0 Implant     | 70-1190-PRA0033 |
|  | Ø3.5/4.3 Implant | 70-1190-PRA0034 |
|  | Ø5.0 Implant     | 70-1190-PRA0035 |
|  | Ø6.0 Implant     | 70-1190-PRA0071 |
|  | Ø7.0 Implant     | 70-1190-PRA0036 |

(UCLA Abutments include Titanium Screw.)

## Stock Abutments

|                                                                                                                         | Description                          | Code                               |
|-------------------------------------------------------------------------------------------------------------------------|--------------------------------------|------------------------------------|
| <b>Glidewell HT™ Implant Titanium Screw and Guide Pin</b>                                                               |                                      |                                    |
| <br>70-1190-PRC0029                    | Ø3.0 Implant                         | 70-1190-PRC0029                    |
| <br>70-1190-PRC0030                    | Ø3.5/4.3 Implant                     | 70-1190-PRC0030                    |
| <br>70-1190-PRC0064<br>70-1190-PRC0031 | Ø5.0/6.0 Implant<br>Ø5.0/7.0 Implant | 70-1190-PRC0064<br>70-1190-PRC0031 |
| <br>70-1190-PRC0047                    | Ø3.0 Implant                         | 70-1190-PRC0047                    |
| <br>70-1190-PRC0048                    | Ø3.5/4.3 Implant                     | 70-1190-PRC0048                    |
| <br>70-1190-PRC0065<br>70-1190-PRC0049 | Ø5.0/6.0 Implant<br>Ø5.0/7.0 Implant | 70-1190-PRC0065<br>70-1190-PRC0049 |

## Multi-Unit Abutments

|                                                                                                        |                            |                 |
|--------------------------------------------------------------------------------------------------------|----------------------------|-----------------|
| <b>Glidewell HT™ Implant Multi-Unit Abutment</b>                                                       |                            |                 |
| <br>70-1190-PRS0011 | 2.5 mmH - Ø3.0 Implant     | 70-1190-PRS0011 |
| <br>70-1190-PRS0014 | 3.5 mmH - Ø3.0 Implant     | 70-1190-PRS0014 |
| <br>70-1190-PRS0009 | 1.5 mmH - Ø3.5/4.3 Implant | 70-1190-PRS0009 |
| <br>70-1190-PRS0012 | 2.5 mmH - Ø3.5/4.3 Implant | 70-1190-PRS0012 |
| <br>70-1190-PRS0015 | 3.5 mmH - Ø3.5/4.3 Implant | 70-1190-PRS0015 |
| <br>70-1190-PRS0010 | 1.5 mmH - Ø5.0 Implant     | 70-1190-PRS0010 |
| <br>70-1190-PRS0013 | 2.5 mmH - Ø5.0 Implant     | 70-1190-PRS0013 |
| <br>70-1190-PRS0016 | 3.5 mmH - Ø5.0 Implant     | 70-1190-PRS0016 |
| <br>70-1190-PRS0017 | 4.5 mmH - Ø5.0 Implant     | 70-1190-PRS0017 |
| <br>70-1190-PRS0019 | 1.5 mmH - Ø6.0 Implant     | 70-1190-PRS0019 |
| <br>70-1190-PRS0020 | 2.5 mmH - Ø6.0 Implant     | 70-1190-PRS0020 |
| <br>70-1190-PRS0021 | 3.5 mmH - Ø6.0 Implant     | 70-1190-PRS0021 |
| <br>70-1190-PRS0022 | 4.5 mmH - Ø6.0 Implant     | 70-1190-PRS0022 |

| Description                                                                                          |                            | Code            |
|------------------------------------------------------------------------------------------------------|----------------------------|-----------------|
| <b>Glidewell HT™ Implant 17° Multi-Unit Abutment</b>                                                 |                            |                 |
| <br>70-1190-PRS0001 | 2.5 mmH - Ø3.5/4.3 Implant | 70-1190-PRS0001 |
| <br>70-1190-PRS0003 | 3.5 mmH - Ø3.5/4.3 Implant | 70-1190-PRS0003 |
| <br>70-1190-PRS0002 | 2.5 mmH - Ø5.0 Implant     | 70-1190-PRS0002 |
| <br>70-1190-PRS0004 | 3.5 mmH - Ø5.0 Implant     | 70-1190-PRS0004 |

|                                                                                                      |                            |                 |
|------------------------------------------------------------------------------------------------------|----------------------------|-----------------|
| <b>Glidewell HT™ Implant 30° Multi-Unit Abutment</b>                                                 |                            |                 |
| <br>70-1190-PRS0005 | 3.5 mmH - Ø3.5/4.3 Implant | 70-1190-PRS0005 |
| <br>70-1190-PRS0007 | 4.5 mmH - Ø3.5/4.3 Implant | 70-1190-PRS0007 |
| <br>70-1190-PRS0006 | 3.5 mmH - Ø5.0 Implant     | 70-1190-PRS0006 |
| <br>70-1190-PRS0008 | 4.5 mmH - Ø5.0 Implant     | 70-1190-PRS0008 |

(Angled Multi-Unit Abutments include Angled Abutment Screw. Prosthetic Screw sold separately; please see page 15.)

|                                                                                                        |                                        |                 |
|--------------------------------------------------------------------------------------------------------|----------------------------------------|-----------------|
| <b>Glidewell HT™ Implant Multi-Unit Coping</b>                                                         |                                        |                 |
| <br>70-1190-PRS0018 | Multi-Unit Coping 4.5mmH, Non-Engaging | 70-1190-PRS0018 |

|                                                                                                        |                  |                 |
|--------------------------------------------------------------------------------------------------------|------------------|-----------------|
| <b>Glidewell HT™ Implant Angled Multi-Unit Abutment Screw</b>                                          |                  |                 |
| <br>70-1190-PRC0032 | Ø3.5/4.3 Implant | 70-1190-PRC0032 |
| <br>70-1190-PRC0033 | Ø5.0 Implant     | 70-1190-PRC0033 |




# Multi-Unit Accessories

The Glidewell HT™ Implant System proudly features Inclusive® and Newport Surgical™ brand multi-unit accessories.

|                                                                                                        | Description                                                    | Code            |
|--------------------------------------------------------------------------------------------------------|----------------------------------------------------------------|-----------------|
| <b>Glidewell HT™ Implant Multi-Unit</b>                                                                |                                                                |                 |
| <br>70-1190-PRC0046   | Titanium Scan Body with Prosthetic Screw                       | 70-1190-PRC0046 |
| <b>Inclusive® Multi-Unit</b>                                                                           |                                                                |                 |
| <br>70-1074-PRS0044   | Temporary Healing Cap with Prosthetic Screw - Unigrip™ Driver  | 70-1074-PRS0044 |
| <br>70-1074-PRS0057   | Titanium Temporary with Prosthetic Screw - Unigrip™ Driver     | 70-1074-PRS0057 |
| <br>70-1072-PRC0043   | Prosthetic Screw - Unigrip™ Driver                             | 70-1072-PRC0043 |
| <br>70-1074-PRS0030   | Gold/Plastic Coping                                            | 70-1074-PRS0030 |
| <br>70-1074-PRS0018  | Gold Bar Coping                                                | 70-1074-PRS0018 |
| <br>70-1074-PRS0069 | Plastic Coping with Prosthetic Screw - Unigrip™ Driver         | 70-1074-PRS0069 |
| <br>70-1074-PRS0022 | Impression Coping, Open-Tray, with Guide Pin - Unigrip™ Driver | 70-1074-PRS0022 |
| <br>70-1074-PRS0021 | Impression Coping, Closed-Tray                                 | 70-1074-PRS0021 |
| <br>70-1074-PRS0016 | Abutment Analog                                                | 70-1074-PRS0016 |
| <br>70-1074-PRS0019 | Guide Pin - Unigrip™ Driver                                    | 70-1074-PRS0019 |
| <i>Unigrip is a trademark of Nobel Biocare Services AG.</i>                                            |                                                                |                 |
| <b>Newport Surgical™ Multi-Unit</b>                                                                    |                                                                |                 |
| <br>70-1071-SRG0234 | Newport Surgical Multi-Unit Driver, 4mm Square                 | 70-1071-SRG0234 |


# LOCATOR® Abutment System

The Glidewell HT™ Implant System proudly features Locator® brand overdenture attachments.
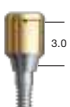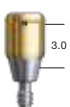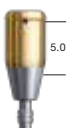
| Description                                                                                            | Code                                                               |
|--------------------------------------------------------------------------------------------------------|--------------------------------------------------------------------|
| <b>LOCATOR® Abutment</b>                                                                               |                                                                    |
| <br>70-1075-ZST0065   | Conical Connection Ø3.0 x 1.0 mmH Cuff [#2212] 70-1075-ZST0065     |
| <br>70-1075-ZST0074   | Conical Connection Ø3.0 x 2.0 mmH Cuff [#2213] 70-1075-ZST0074     |
| <br>70-1075-ZST0066   | Conical Connection Ø3.0 x 3.0 mmH Cuff [#2214] 70-1075-ZST0066     |
| <br>70-1075-ZST0032   | Conical Connection Ø3.5/4.3 x 1.0 mmH Cuff [#2069] 70-1075-ZST0032 |
| <br>70-1075-ZST0034  | Conical Connection Ø3.5/4.3 x 3.0 mmH Cuff [#2071] 70-1075-ZST0034 |
| <br>70-1075-ZST0036 | Conical Connection Ø3.5/4.3 x 5.0 mmH Cuff [#2073] 70-1075-ZST0036 |
| <br>70-1075-ZST0071 | Conical Connection Ø5.0 x 1.0 mmH Cuff [#2075] 70-1075-ZST0071     |
| <br>70-1075-ZST0038 | Conical Connection Ø5.0 x 3.0 mmH Cuff [#2077] 70-1075-ZST0038     |
| <br>70-1075-ZST0040 | Conical Connection Ø5.0 x 5.0 mmH Cuff [#2079] 70-1075-ZST0040     |

# LOCATOR Abutment System

The Glidewell HT Implant System proudly features Locator brand overdenture attachments.
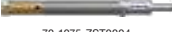
| Description                                                                                                | Code                                                                                                                                                                                                                                                                                |
|------------------------------------------------------------------------------------------------------------|-------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|
| <b>LOCATOR®</b>                                                                                            |                                                                                                                                                                                                                                                                                     |
|  <p>70-1075-ZST0024</p>   | <b>Male Processing Package (2-pack) [#8519-2]</b> 70-1075-ZST0024<br><i>(Includes: Denture Cap with Black Processing Male, Clear Nylon Male (5 lbs), Pink (LR) Less-Retentive Nylon Male (3 lbs), Blue Extra-Light-Retention Nylon Male (1.5 lbs), and White Block-out Spacer.)</i> |
|  <p>70-1075-ZST0025</p>   | <b>Extended Range Male Processing Package (2-pack) [#8540-2]</b> 70-1075-ZST0025<br><i>(Includes: Denture Cap with Black Processing Male, Extended-Range Replacement Males: Green (3-4 lbs), Orange (2 lbs), Red (0.5-1.5 lbs), and White Block-out Spacer.)</i>                    |
|  <p>70-1075-ZST0007</p>   | <b>Impression Coping (4-pack) [#8505]</b> 70-1075-ZST0007                                                                                                                                                                                                                           |
|  <p>70-1075-ZST0009</p>   | <b>Female Analog 4 mmD (4-pack) [#8530]</b> 70-1075-ZST0009                                                                                                                                                                                                                         |
|  <p>70-1075-ZST0005</p>   | <b>Core Tool Abutment Holder Sleeve (4-pack) [#8394]</b> 70-1075-ZST0005                                                                                                                                                                                                            |
|  <p>70-1075-ZST0026</p>   | <b>Angle Measurement Guide [#9530]</b> 70-1075-ZST0026                                                                                                                                                                                                                              |
|  <p>70-1075-ZST0010</p>   | <b>Processing Spacer (4-pack) [#8569]</b> 70-1075-ZST0010                                                                                                                                                                                                                           |
|  <p>70-1075-ZST0004</p> | <b>Core Tool (Complete) [#8393]</b> 70-1075-ZST0004                                                                                                                                                                                                                                 |
|  <p>70-1075-ZST0006</p> | <b>Core Tool Male Removal (Tip End) [#8397]</b> 70-1075-ZST0006                                                                                                                                                                                                                     |
|  <p>70-1075-ZST0008</p> | <b>Parallel Post (4-pack) [#8517]</b> 70-1075-ZST0008                                                                                                                                                                                                                               |
|  <p>70-1075-ZST0003</p> | <b>Core Tool Abutment Driver (Gold End) [#8390]</b> 70-1075-ZST0003                                                                                                                                                                                                                 |
|  <p>70-1075-ZST0001</p> | <b>35 Ncm Torque Wrench Insert Driver (15 mm) [#8260]</b> 70-1075-ZST0001                                                                                                                                                                                                           |
|  <p>70-1075-ZST0002</p> | <b>35 Ncm Torque Wrench Insert Driver (21 mm) [#8280]</b> 70-1075-ZST0002                                                                                                                                                                                                           |
|  <p>70-1075-ZST0020</p> | <b>Latch Type Torque Wrench Insert Driver (23 mm) [#8913]</b> 70-1075-ZST0020                                                                                                                                                                                                       |
|  <p>70-1075-ZST0021</p> | <b>Latch Type Torque Wrench Insert Driver (29 mm) [#8914]</b> 70-1075-ZST0021                                                                                                                                                                                                       |
|  <p>70-1075-ZST0022</p> | <b>Square Drive Torque Wrench Insert Driver (15 mm) [#8926]</b> 70-1075-ZST0022                                                                                                                                                                                                     |
|  <p>70-1075-ZST0023</p> | <b>Square Drive Torque Wrench Insert Driver (21 mm) [#8927]</b> 70-1075-ZST0023                                                                                                                                                                                                     |

LOCATOR is a registered trademark of Zest Anchors, LLC.

# Standard Surgical Instrumentation Kits

## Features and Benefits

The instrumentation kits allow the clinician to easily organize, store and transport the tooling of the Glidewell HT™ Implant System. The autoclavable kits have been designed for ease of use, with component labels for identification.

The Glidewell HT Implant Surgical Kit includes instrumentation that is machined from corrosion-resistant, surgical stainless steel, and features standard connectivity. The color-coded surgical kit clearly maps the drilling protocol for each implant. Drills are arranged from left to right in order of increasing diameter, with color-coded bands on the shaping drills to indicate the corresponding diameter of the Glidewell HT Implant. For added safety, each shaping drill is implant-specific to precisely determine both diameter and depth of the osteotomy. **Please note that Ø6.0 mm and Ø7.0 mm Shaping Drills, Screw Tap, Direction Indicator, and Driver are sold separately (refer to diagram, page 19).**


*All instrumentation is manufactured in the U.S.A., Switzerland, Denmark or Germany. For specific country of origin, please refer to the individual product label.*


## Sterility

Please refer to the Glidewell HT Implant System Surgical Manual for cleaning and sterilization instructions.

| Description                        | Code            |
|------------------------------------|-----------------|
| Glidewell HT™ Implant Surgical Kit | 70-1071-SRG0396 |



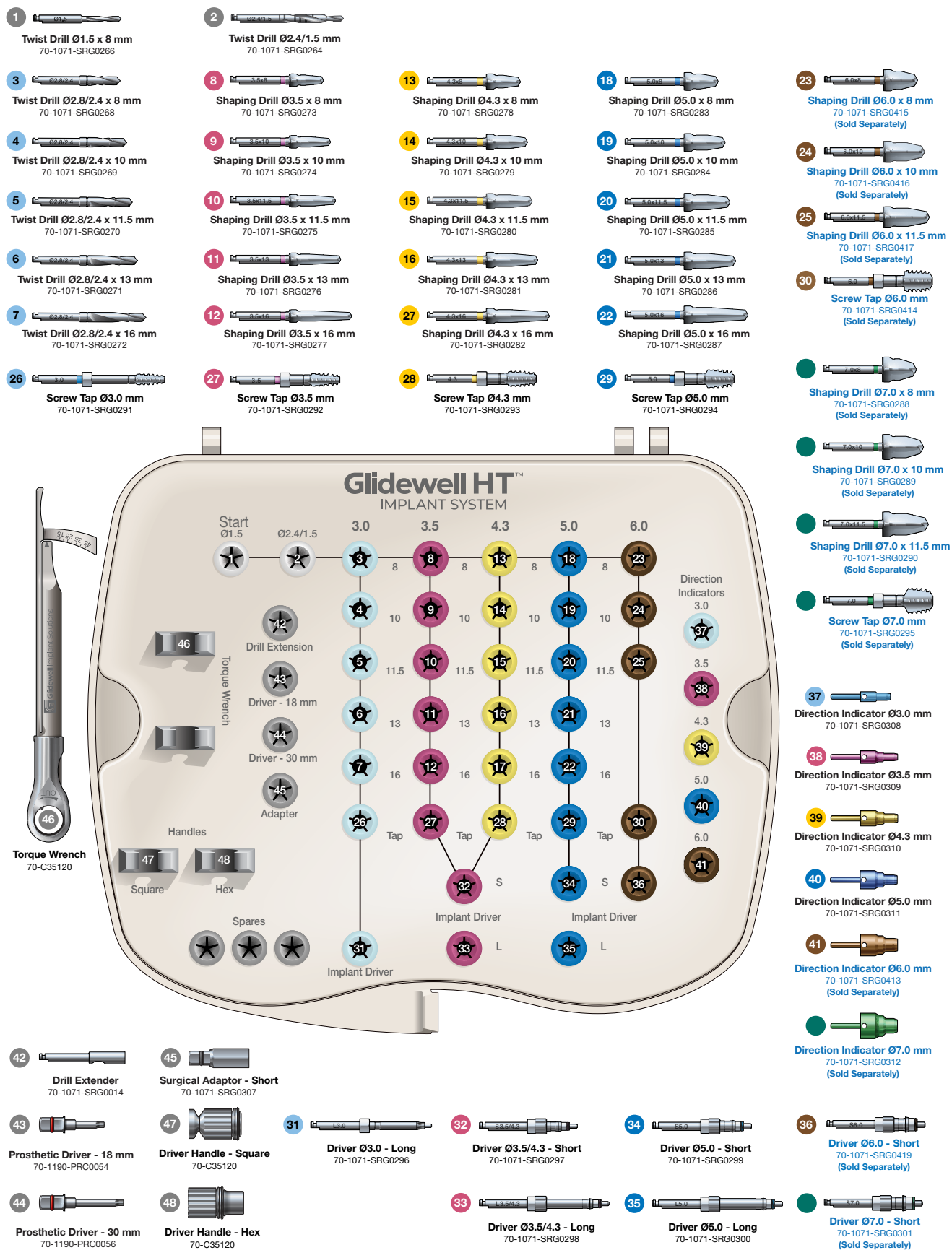

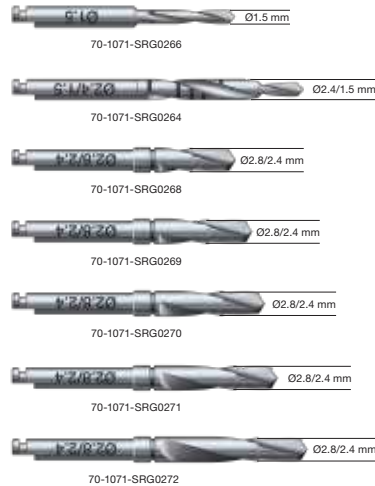

## Glidewell HT™ Implant Twist Drill

| Description        | Code            |
|--------------------|-----------------|
| Ø1.5 x 8 mm        | 70-1071-SRG0266 |
| Ø2.4/1.5 mm        | 70-1071-SRG0264 |
| Ø2.8/2.4 x 8 mm    | 70-1071-SRG0268 |
| Ø2.8/2.4 x 10 mm   | 70-1071-SRG0269 |
| Ø2.8/2.4 x 11.5 mm | 70-1071-SRG0270 |
| Ø2.8/2.4 x 13 mm   | 70-1071-SRG0271 |
| Ø2.8/2.4 x 16 mm   | 70-1071-SRG0272 |



## Glidewell HT™ Implant Shaping Drill

|                |                 |
|----------------|-----------------|
| Ø3.5 x 8 mm    | 70-1071-SRG0273 |
| Ø3.5 x 10 mm   | 70-1071-SRG0274 |
| Ø3.5 x 11.5 mm | 70-1071-SRG0275 |
| Ø3.5 x 13 mm   | 70-1071-SRG0276 |
| Ø3.5 x 16 mm   | 70-1071-SRG0277 |
| Ø4.3 x 8 mm    | 70-1071-SRG0278 |
| Ø4.3 x 10 mm   | 70-1071-SRG0279 |
| Ø4.3 x 11.5 mm | 70-1071-SRG0280 |
| Ø4.3 x 13 mm   | 70-1071-SRG0281 |
| Ø4.3 x 16 mm   | 70-1071-SRG0282 |
| Ø5.0 x 8 mm    | 70-1071-SRG0283 |
| Ø5.0 x 10 mm   | 70-1071-SRG0284 |



70-1071-SRG0285



70-1071-SRG0286



70-1071-SRG0287



70-1071-SRG0415



70-1071-SRG0416



70-1071-SRG0417



70-1071-SRG0288



70-1071-SRG0289



70-1071-SRG0290

## Glidewell HT™ Implant Shaping Drill

| Description    | Code            |
|----------------|-----------------|
| Ø5.0 x 11.5 mm | 70-1071-SRG0285 |
| Ø5.0 x 13 mm   | 70-1071-SRG0286 |
| Ø5.0 x 16 mm   | 70-1071-SRG0287 |
| Ø6.0 x 8 mm    | 70-1071-SRG0415 |
| Ø6.0 x 10 mm   | 70-1071-SRG0416 |
| Ø6.0 x 11.5 mm | 70-1071-SRG0417 |
| Ø7.0 x 8 mm    | 70-1071-SRG0288 |
| Ø7.0 x 10 mm   | 70-1071-SRG0289 |
| Ø7.0 x 11.5 mm | 70-1071-SRG0290 |



70-1071-SRG0291



70-1071-SRG0292



70-1071-SRG0293



70-1071-SRG0294




70-1071-SRG0414



70-1071-SRG0295


## Glidewell HT™ Implant Screw Tap

| Description | Code            |
|-------------|-----------------|
| Ø3.0 mm     | 70-1071-SRG0291 |
| Ø3.5 mm     | 70-1071-SRG0292 |
| Ø4.3 mm     | 70-1071-SRG0293 |
| Ø5.0 mm     | 70-1071-SRG0294 |
| Ø6.0 mm     | 70-1071-SRG0414 |
| Ø7.0 mm     | 70-1071-SRG0295 |

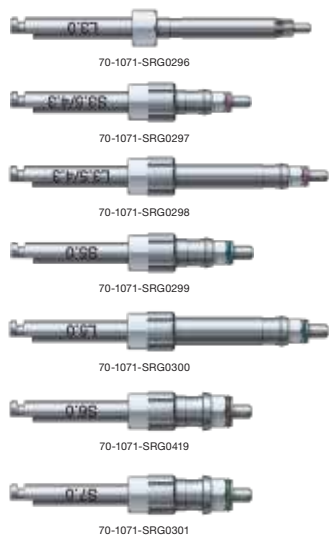


## Glidewell HT™ Implant Driver

| Description    | Code            |
|----------------|-----------------|
| Ø3.0 Long      | 70-1071-SRG0296 |
| Ø3.5/4.3 Short | 70-1071-SRG0297 |
| Ø3.5/4.3 Long  | 70-1071-SRG0298 |
| Ø5.0 Short     | 70-1071-SRG0299 |
| Ø5.0 Long      | 70-1071-SRG0300 |
| Ø6.0 Short     | 70-1071-SRG0419 |
| Ø7.0 Short     | 70-1071-SRG0301 |

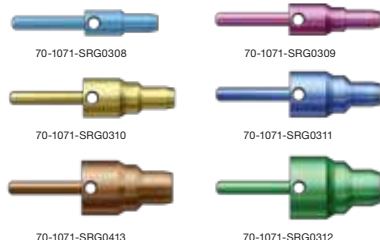


## Glidewell HT™ Implant Direction Indicator

|         |                 |
|---------|-----------------|
| Ø3.0 mm | 70-1071-SRG0308 |
| Ø3.5 mm | 70-1071-SRG0309 |
| Ø4.3 mm | 70-1071-SRG0310 |
| Ø5.0 mm | 70-1071-SRG0311 |
| Ø6.0 mm | 70-1071-SRG0413 |
| Ø7.0 mm | 70-1071-SRG0312 |

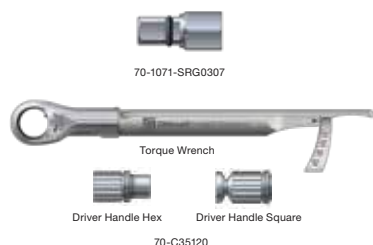

|                                                |                 |
|------------------------------------------------|-----------------|
| Glidewell HT™ Implant Surgical Adaptor - Short | 70-1071-SRG0307 |
|------------------------------------------------|-----------------|

|                                                                |           |
|----------------------------------------------------------------|-----------|
| Glidewell HT™ Torque Wrench with Driver Handle, Hex and Square | 70-C35120 |
|----------------------------------------------------------------|-----------|

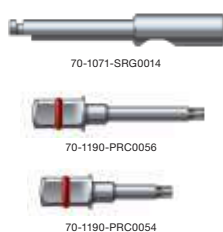

|                           |                 |
|---------------------------|-----------------|
| Inclusive® Drill Extender | 70-1071-SRG0014 |
|---------------------------|-----------------|

|                                         |                 |
|-----------------------------------------|-----------------|
| Glidewell HT™ Prosthetic Driver - 30 mm | 70-1190-PRC0056 |
|-----------------------------------------|-----------------|

|                                         |                 |
|-----------------------------------------|-----------------|
| Glidewell HT™ Prosthetic Driver - 18 mm | 70-1190-PRC0054 |
|-----------------------------------------|-----------------|



|                                             |                 |
|---------------------------------------------|-----------------|
| Glidewell HT™ Implant Radiographic Template | 70-1071-RAD0001 |
|---------------------------------------------|-----------------|

# Guided Surgical Instrumentation Kit

## Features and Benefits

The Glidewell HT™ Implant Guided Surgery System allows clinicians to deliver Glidewell HT Implants with the highest level of predictability, accuracy, and efficiency. Designed for ease of use, including enhancements based on eight years of clinical feedback, this system works in conjunction with digital treatment planning and case-specific surgical guides.

The Glidewell HT™ Implant Guided Surgical Kit includes instrumentation that is machined from high-quality, corrosion-resistant, surgical stainless steel. Color-coded bands correspond to implant diameter, matching the color of the implant carrier. Tissue punches match the diameter of the prescribed implant. Drills are diameter-specific and feature a flange stop for depth control. For added safety, each shaping drill is implant-specific to precisely determine both diameter and depth of the osteotomy. Drills are arranged from left to right in order of increasing diameter, following the recommended drilling sequence.

*All instrumentation is manufactured in the U.S.A., Switzerland, Denmark or Israel. For specific country of origin, please refer to the individual product label.*



## Sterility

Please refer to the Glidewell HT Implant System Surgical Manual for cleaning and sterilization instructions.


| Description                               | Code            |
|-------------------------------------------|-----------------|
| Glidewell HT™ Implant Guided Surgical Kit | 70-1071-SRG0395 |



# Guided Surgical Kit Overview

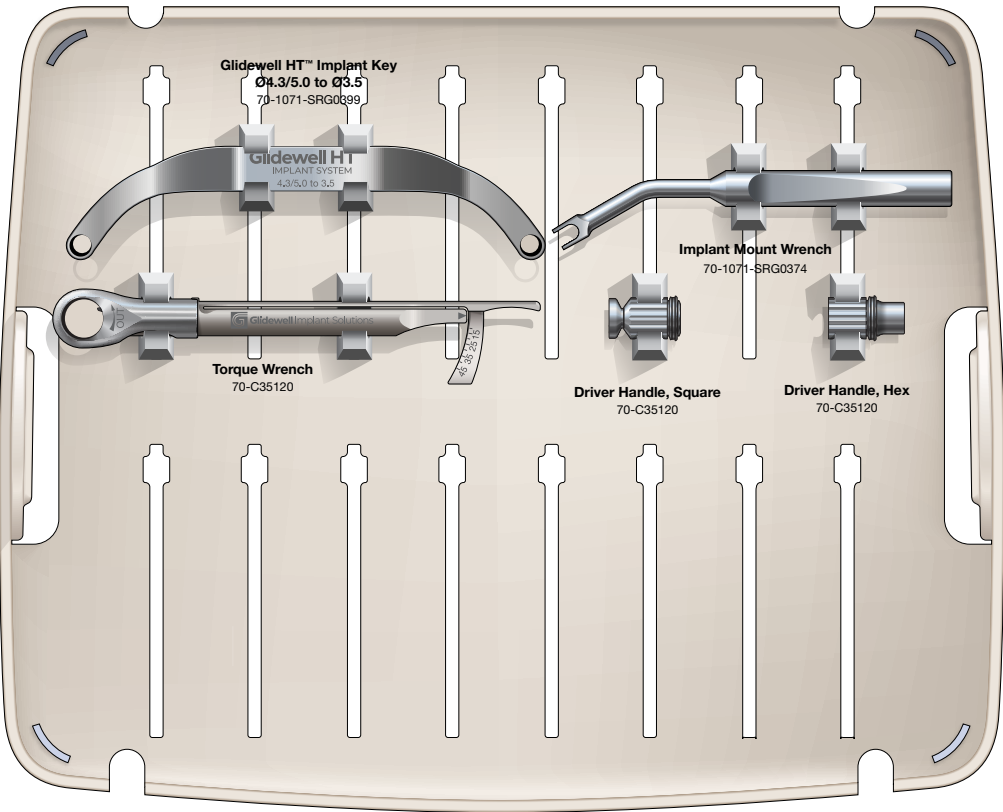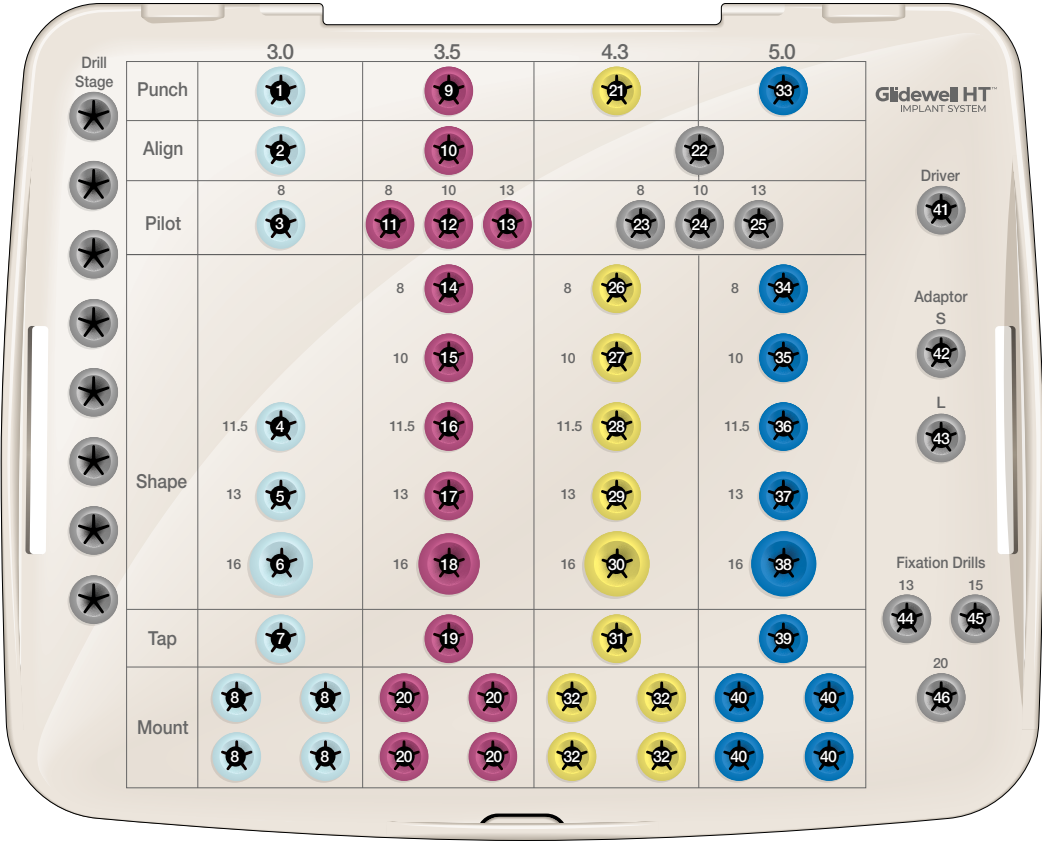

# Guided Surgical Kit Overview

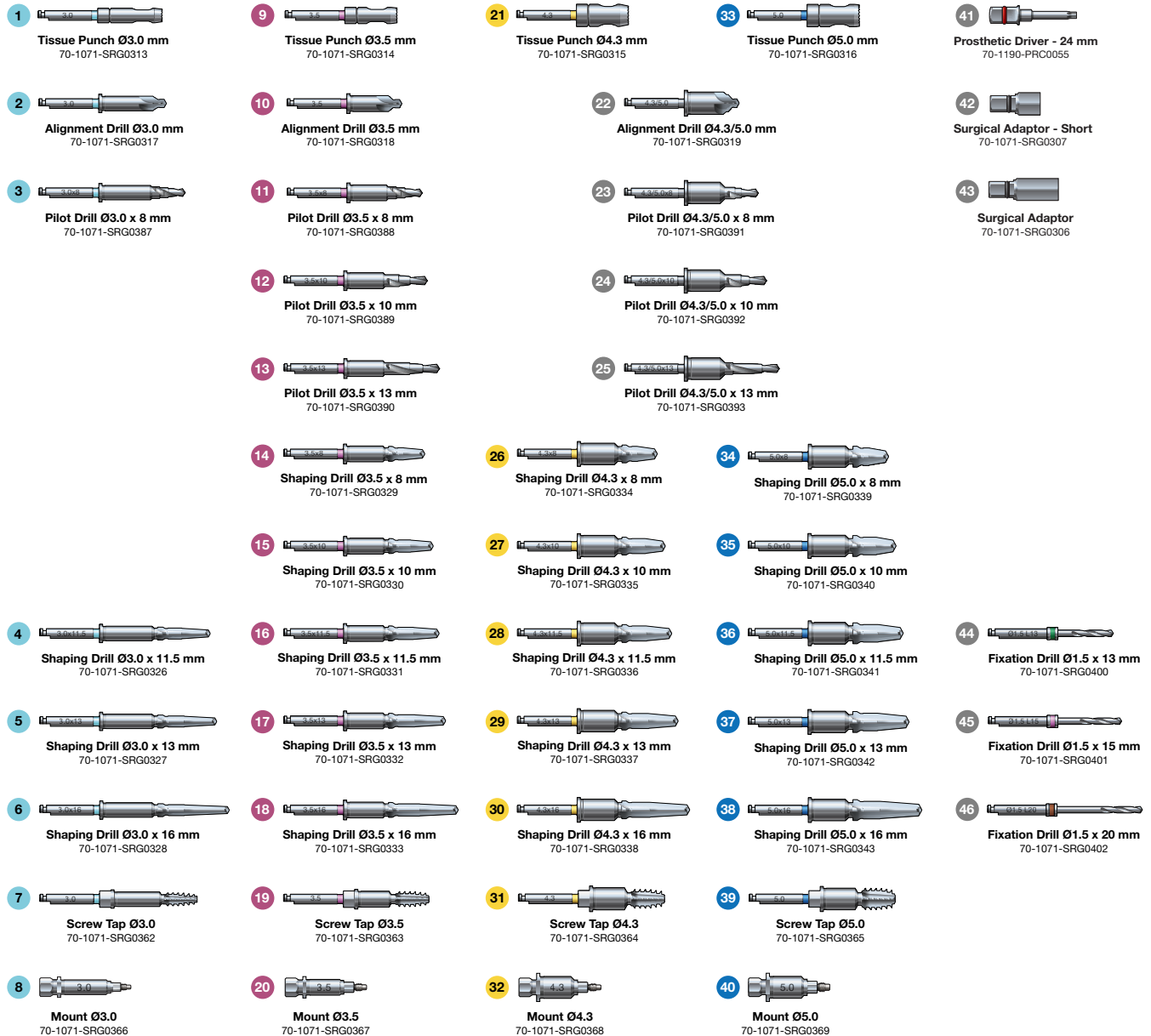



70-1071-SRG0313



70-1071-SRG0314



70-1071-SRG0315



70-1071-SRG0316

Description

Code

## Glidewell HT™ Implant Guided Tissue Punch

|         |                 |
|---------|-----------------|
| Ø3.0 mm | 70-1071-SRG0313 |
| Ø3.5 mm | 70-1071-SRG0314 |
| Ø4.3 mm | 70-1071-SRG0315 |
| Ø5.0 mm | 70-1071-SRG0316 |



70-1071-SRG0317



70-1071-SRG0318



70-1071-SRG0319

## Glidewell HT™ Implant Guided Alignment Drill

|             |                 |
|-------------|-----------------|
| Ø3.0 mm     | 70-1071-SRG0317 |
| Ø3.5 mm     | 70-1071-SRG0318 |
| Ø4.3/5.0 mm | 70-1071-SRG0319 |



70-1071-SRG0387



70-1071-SRG0388



70-1071-SRG0389



70-1071-SRG0390



70-1071-SRG0391




70-1071-SRG0392



70-1071-SRG0393

## Glidewell HT™ Implant Guided Pilot Drill

|                  |                 |
|------------------|-----------------|
| Ø3.0 x 8 mm      | 70-1071-SRG0387 |
| Ø3.5 x 8 mm      | 70-1071-SRG0388 |
| Ø3.5 x 10 mm     | 70-1071-SRG0389 |
| Ø3.5 x 13 mm     | 70-1071-SRG0390 |
| Ø4.3/5.0 x 8 mm  | 70-1071-SRG0391 |
| Ø4.3/5.0 x 10 mm | 70-1071-SRG0392 |
| Ø4.3/5.0 x 13 mm | 70-1071-SRG0393 |
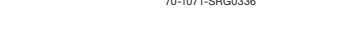
|                                                                                                        | Description    | Code            |
|--------------------------------------------------------------------------------------------------------|----------------|-----------------|
| <b>Glidewell HT™ Implant Guided Shaping Drill</b>                                                      |                |                 |
| <br>70-1071-SRG0326   | Ø3.0 x 11.5 mm | 70-1071-SRG0326 |
| <br>70-1071-SRG0327   | Ø3.0 x 13 mm   | 70-1071-SRG0327 |
| <br>70-1071-SRG0328   | Ø3.0 x 16 mm   | 70-1071-SRG0328 |
| <br>70-1071-SRG0329   | Ø3.5 x 8 mm    | 70-1071-SRG0329 |
| <br>70-1071-SRG0330   | Ø3.5 x 10 mm   | 70-1071-SRG0330 |
| <br>70-1071-SRG0331   | Ø3.5 x 11.5 mm | 70-1071-SRG0331 |
| <br>70-1071-SRG0332  | Ø3.5 x 13 mm   | 70-1071-SRG0332 |
| <br>70-1071-SRG0333 | Ø3.5 x 16 mm   | 70-1071-SRG0333 |
| <br>70-1071-SRG0334 | Ø4.3 x 8 mm    | 70-1071-SRG0334 |
| <br>70-1071-SRG0335 | Ø4.3 x 10 mm   | 70-1071-SRG0335 |
| <br>70-1071-SRG0336 | Ø4.3 x 11.5 mm | 70-1071-SRG0336 |
| <br>70-1071-SRG0337 | Ø4.3 x 13 mm   | 70-1071-SRG0337 |
| <br>70-1071-SRG0338 | Ø4.3 x 16 mm   | 70-1071-SRG0338 |

## Guided Drills



70-1071-SRG0339



70-1071-SRG0340



70-1071-SRG0341



70-1071-SRG0342



70-1071-SRG0343

Description

Code

### Glidewell HT™ Implant Guided Shaping Drills

Ø5.0 x 8 mm

70-1071-SRG0339

Ø5.0 x 10 mm

70-1071-SRG0340

Ø5.0 x 11.5 mm

70-1071-SRG0341

Ø5.0 x 13 mm

70-1071-SRG0342

Ø5.0 x 16 mm

70-1071-SRG0343

## Guided Taps & Accessories



70-1071-SRG0362



70-1071-SRG0363



70-1071-SRG0364



70-1071-SRG0365

### Glidewell HT™ Implant Guided Screw Tap

Ø3.0 mm

70-1071-SRG0362

Ø3.5 mm

70-1071-SRG0363

Ø4.3 mm

70-1071-SRG0364

Ø5.0 mm

70-1071-SRG0365

Description

Code

## Glidewell HT™ Implant Guided Mount



70-1071-SRG0366

Ø3.0 Implant

70-1071-SRG0366



70-1071-SRG0367

Ø3.5 Implant

70-1071-SRG0367



70-1071-SRG0368

Ø4.3 Implant

70-1071-SRG0368



70-1071-SRG0369

Ø5.0 Implant

70-1071-SRG0369

(Guided Mounts include captured screw.)



70-1071-SRG0307

Glidewell HT™ Implant Surgical Adaptor - Short

70-1071-SRG0307



70-1071-SRG0306

Glidewell HT™ Implant Surgical Adaptor

70-1071-SRG0306



Torque Wrench



Driver Handle Hex



Driver Handle Square

70-C35120

Glidewell HT™ Torque Wrench with Driver Handle, Hex and Square 70-C35120



70-1071-SRG0374

Glidewell HT™ Implant Mount Wrench

70-1071-SRG0374



70-1190-PRC0055

Glidewell HT™ Prosthetic Driver - 24 mm

70-1190-PRC0055



70-1071-SRG0399

Glidewell HT™ Implant Key - Ø4.3/5.0 to Ø3.5

70-1071-SRG0399

# Prosthetic Instrumentation Kit

## Features and Benefits

The prosthetic kit is designed specifically for restorative doctors, organizing all of the necessary prosthetic tooling to help simplify restorative procedures and instrument sterilization.

*All instrumentation is manufactured in the U.S.A., Switzerland or Denmark. For specific country of origin, please refer to the individual product label.*

## Sterility

Please refer to the Glidewell HT Implant System Restorative Manual for cleaning and sterilization instructions.

| Description               | Code            |
|---------------------------|-----------------|
| Glidewell™ Prosthetic Kit | 70-1071-PRS0004 |



# Prosthetic Instrumentation Kit Overview

**Glidewell HT™**  
IMPLANT SYSTEM

**11**   
**Multi-Unit Driver**  
70-1071-SRG0234


**Glidewell IH™**  
IMPLANT SYSTEM

**1**   
**Prosthetic Driver - 18 mm**  
70-1190-PRC0054

**2**   
**Prosthetic Driver - 24 mm**  
70-1190-PRC0055

**3**   
**Prosthetic Driver - 30 mm**  
70-1190-PRC0056

**4**   
**Prosthetic Driver - 38 mm**  
70-1190-PRC0057

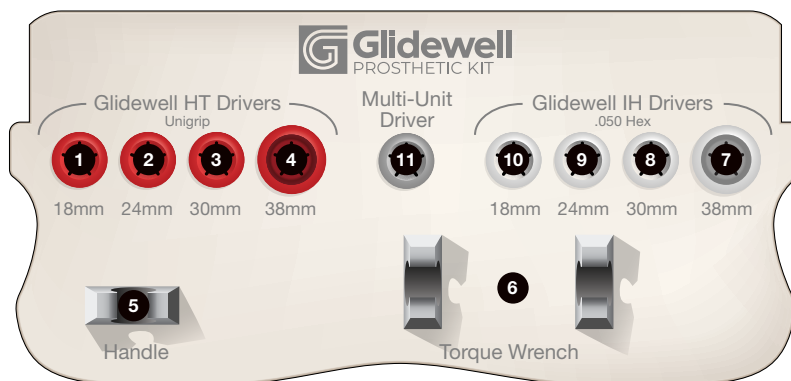


**10**   
**Prosthetic Driver, .050 Hex - 18 mm**  
70-1195-PRC0001

**9**   
**Prosthetic Driver, .050 Hex - 24 mm**  
70-1195-PRC0002

**8**   
**Prosthetic Driver, .050 Hex - 30 mm**  
70-1195-PRC0003

**7**   
**Prosthetic Driver, .050 Hex - 38 mm**  
70-1195-PRC0004

**5**   
**Driver Handle, Square**  
70-C35119

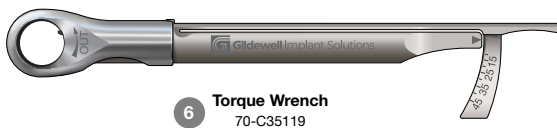**6**   
**Torque Wrench**  
70-C35119

Description

Code



Torque Wrench



Driver Handle Square

70-C35119



70-1071-SRG0234



70-1190-PRC0054



70-1190-PRC0055



70-1190-PRC0056



70-1190-PRC0057



70-1195-PRC0001



70-1195-PRC0002



70-1195-PRC0003



70-1195-PRC0004

Glidewell HT™ Torque Wrench with Driver Handle, Square

70-C35119

Newport Surgical™ Multi-Unit Driver, 4 mm Square

70-1071-SRG0234

Glidewell HT™ Prosthetic Driver - 18 mm

70-1190-PRC0054

Glidewell HT™ Prosthetic Driver - 24 mm

70-1190-PRC0055

Glidewell HT™ Prosthetic Driver - 30 mm

70-1190-PRC0056

Glidewell HT™ Prosthetic Driver - 38 mm

70-1190-PRC0057

Glidewell IH™ Prosthetic Driver - 18 mm

70-1195-PRC0001

Glidewell IH™ Prosthetic Driver - 24 mm

70-1195-PRC0002

Glidewell IH™ Prosthetic Driver - 30 mm

70-1195-PRC0003

Glidewell IH™ Prosthetic Driver - 38 mm



70-1195-PRC0004

## Product Packaging

Glidewell HT™ Implants are shipped sterile. They should not be resterilized. They are for single use only, prior to the expiration date. Do not use implants if the packaging has been compromised or previously opened. Do not handle implant surfaces directly. Users are advised to visually inspect packaging to ensure seals and contents are intact prior to use. Please refer to the individual product label for all relevant product information and cautions.

### Explanation of Label Codes:

1. Official product description
2. Date of Manufacture (YYYY-MM-DD)
3. Catalog Number
4. Lot/Batch Number
5. By Prescription Only
6. Do Not Use if Package is Damaged
7. Sterile with Gamma Radiation
8. Use-by Date
9. Consult Instructions for Use
10. Do not Re-use
11. Do not Resterilize
12. Manufacturer
13. Country of origin
14. Unique Device Identification (UDI)
15. Non-Sterile
16. QR code for IFU website

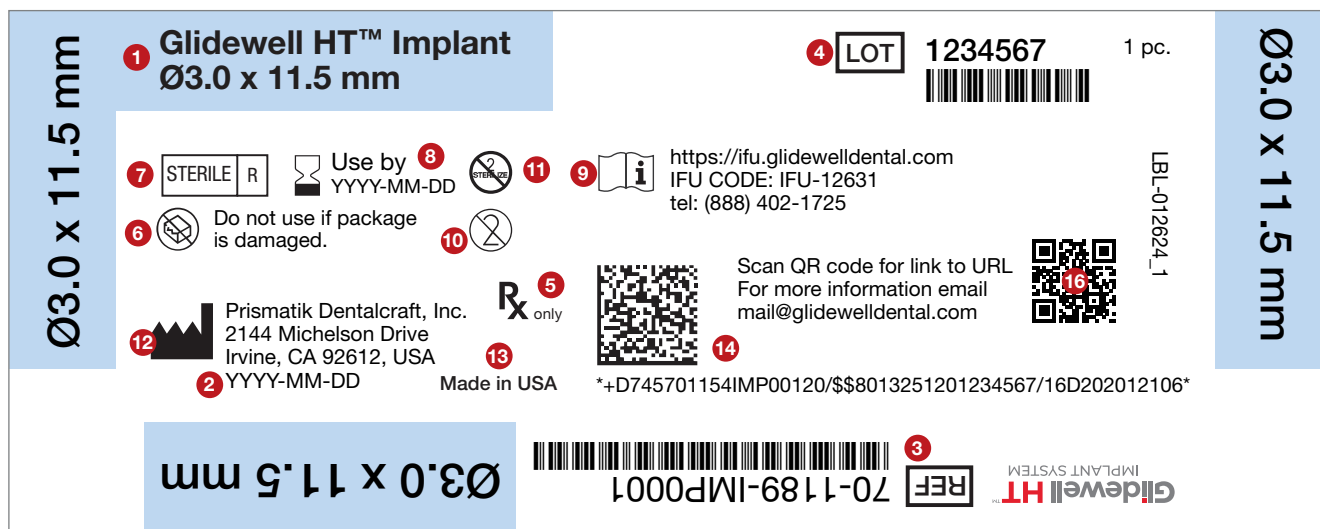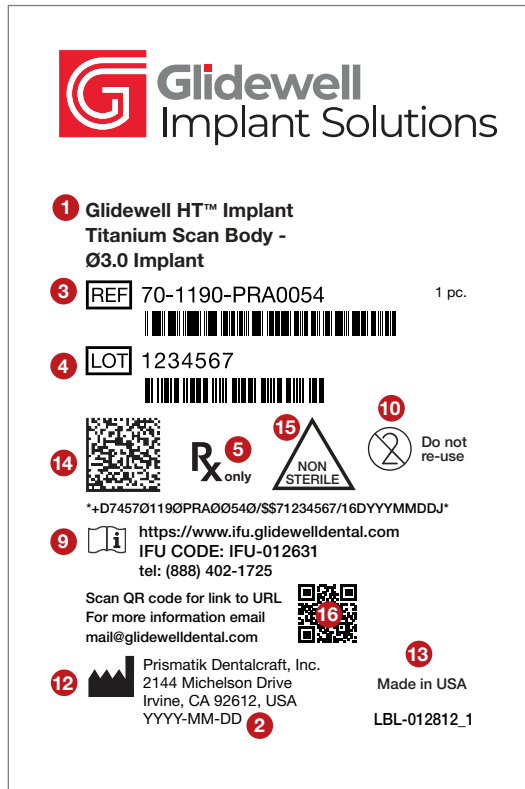

### Ordering Information

In the U.S.: Order at [glidewelldirect.com](http://glidewelldirect.com) or call Glidewell Direct at 888-303-3975. Our product specialists are committed to answering questions in a timely fashion to ensure your ordering is easy and efficient. We are available Monday—Friday from 6 a.m.—5 p.m. (*Pacific*).

### Shipping Policy

- Orders placed after 3 p.m. (*Pacific*) will be processed on the following business day. Business days do not include Saturdays, Sundays, or U.S. holidays.
- Online shopping cart available to U.S. customers only.

### Terms

All accounts are payable within 30 days of invoice date. Accounts not paid within the stated terms will be subject to COD status and a late charge of 2 percent of the unpaid balance. We accept American Express, Visa, MasterCard, and Discover. All prices are subject to change without notice.

### Product Return Policy

Products may be returned at the customer's expense for credit within 30 days of invoice date. All returned products must meet the following conditions:

- A copy of the original invoice must accompany the products.
- Products must be packaged to arrive at the seller's facility undamaged.
- Discontinued, obsolete, expired, damaged, or opened items will not be accepted for return.
- Amount credited will be based on invoice price, less 15 percent for restocking fee.
- Shipping charges are the responsibility of the customer and will not be credited.

### Product & Pricing Changes

Because products and equipment are continually undergoing refinement in design and manufacturing methods, we reserve the right to improve, modify, or discontinue products and equipment or change pricing at any time without incurring any obligation and without prior notice.

### Warranty

Limited Warranty—Prismatik Dentalcraft, Inc.

Prismatik Dentalcraft, Inc. ("Prismatik"), is the manufacturer of dental products (the "product"), including Glidewell HT™ Implants ("implants"). Prismatik warrants the Glidewell HT Implant for the life of the patient originally receiving the implant from the date of placement, and for a period of six (6) months for ceramic blanks and any other product ("the warranty period"). Prismatik will at its option replace or refund the purchase price of any product, to the original purchaser ("user"), that is returned due to defects in material and manufacture.

NO GUARANTEE OR WARRANTY IS IMPLIED OTHER THAN EXPRESSLY STATED, INCLUDING ANY IMPLIED WARRANTY OF MERCHANTABILITY OR FITNESS FOR A PARTICULAR PURPOSE.

Prismatik shall not be liable for any incidental or consequential damages, whether foreseeable or not, caused by defects in the product or dental devices produced using said product. User is responsible for determining the suitability of the product for user's application. If this product is defective within the warranty period, user's exclusive remedy and Prismatik's sole obligation shall be replacement or refund of the purchase price of the product. For replacement or refund under this warranty, the original purchaser shall send the product at its own expense, postage prepaid, to the seller.

| CODE                                       | DESCRIPTION                       | PAGE | CODE                                               | DESCRIPTION                                           | PAGE   |
|--------------------------------------------|-----------------------------------|------|----------------------------------------------------|-------------------------------------------------------|--------|
| <b>IMPLANTS</b>                            |                                   |      |                                                    |                                                       |        |
| 70-1189-IMP0001                            | .....Implant Ø3.0 x 11.5 mm       | 6    | 70-1071-SRG0415                                    | .....Shaping Drill Ø6.0 x 8 mm                        | 21     |
| 70-1189-IMP0002                            | .....Implant Ø3.0 x 13 mm         | 6    | 70-1071-SRG0416                                    | .....Shaping Drill Ø6.0 x 10 mm                       | 21     |
| 70-1189-IMP0003                            | .....Implant Ø3.0 x 16 mm         | 6    | 70-1071-SRG0417                                    | .....Shaping Drill Ø6.0 x 11.5 mm                     | 21     |
| 70-1189-IMP0004                            | .....Implant Ø3.5 x 8 mm          | 6    | 70-1071-SRG0288                                    | .....Shaping Drill Ø7.0 x 8 mm                        | 21     |
| 70-1189-IMP0005                            | .....Implant Ø3.5 x 10 mm         | 6    | 70-1071-SRG0289                                    | .....Shaping Drill Ø7.0 x 10 mm                       | 21     |
| 70-1189-IMP0006                            | .....Implant Ø3.5 x 11.5 mm       | 6    | 70-1071-SRG0290                                    | .....Shaping Drill Ø7.0 x 11.5 mm                     | 21     |
| 70-1189-IMP0007                            | .....Implant Ø3.5 x 13 mm         | 6    | 70-1071-SRG0264                                    | .....Twist Drill Ø2.4/1.5 mm                          | 20     |
| 70-1189-IMP0008                            | .....Implant Ø3.5 x 16 mm         | 6    | 70-1071-SRG0266                                    | .....Twist Drill Ø1.5 mm x 8 mm                       | 20     |
| 70-1189-IMP0009                            | .....Implant Ø4.3 x 8 mm          | 6    | 70-1071-SRG0268                                    | .....Twist Drill Ø2.8/2.4 mm x 8 mm                   | 20     |
| 70-1189-IMP0010                            | .....Implant Ø4.3 x 10 mm         | 6    | 70-1071-SRG0269                                    | .....Twist Drill Ø2.8/2.4 mm x 10 mm                  | 20     |
| 70-1189-IMP0011                            | .....Implant Ø4.3 x 11.5 mm       | 6    | 70-1071-SRG0270                                    | .....Twist Drill Ø2.8/2.4 mm x 11.5 mm                | 20     |
| 70-1189-IMP0012                            | .....Implant Ø4.3 x 13 mm         | 6    | 70-1071-SRG0271                                    | .....Twist Drill Ø2.8/2.4 mm x 13 mm                  | 20     |
| 70-1189-IMP0013                            | .....Implant Ø4.3 x 16 mm         | 6    | 70-1071-SRG0272                                    | .....Twist Drill Ø2.8/2.4 mm x 16 mm                  | 20     |
| 70-1189-IMP0014                            | .....Implant Ø5.0 x 8 mm          | 6    | <b>INSTRUMENTATION - DRIVERS &amp; ACCESSORIES</b> |                                                       |        |
| 70-1189-IMP0015                            | .....Implant Ø5.0 x 10 mm         | 6    | 70-C35119                                          | .....Torque Wrench with Driver Handle, Square         | 31     |
| 70-1189-IMP0016                            | .....Implant Ø5.0 x 11.5 mm       | 6    | 70-C35120                                          | .....Torque Wrench with Driver Handle, Hex and Square | 22, 29 |
| 70-1189-IMP0017                            | .....Implant Ø5.0 x 13 mm         | 6    | 70-1071-SRG0014                                    | .....Drill Extender                                   | 22     |
| 70-1189-IMP0018                            | .....Implant Ø5.0 x 16 mm         | 6    | 70-1071-SRG0296                                    | .....Driver Ø3.0, Long                                | 22     |
| 70-1189-IMP0022                            | .....Implant Ø6.0 x 8 mm          | 6    | 70-1071-SRG0297                                    | .....Driver Ø3.5/4.3, Short                           | 22     |
| 70-1189-IMP0023                            | .....Implant Ø6.0 x 10 mm         | 6    | 70-1071-SRG0298                                    | .....Driver Ø3.5/4.3, Long                            | 22     |
| 70-1189-IMP0024                            | .....Implant Ø6.0 x 11.5 mm       | 6    | 70-1071-SRG0299                                    | .....Driver Ø5.0, Short                               | 22     |
| 70-1189-IMP0019                            | .....Implant Ø7.0 x 8 mm          | 6    | 70-1071-SRG0300                                    | .....Driver Ø5.0, Long                                | 22     |
| 70-1189-IMP0020                            | .....Implant Ø7.0 x 10 mm         | 6    | 70-1071-SRG0419                                    | .....Driver Ø6.0, Short                               | 22     |
| 70-1189-IMP0021                            | .....Implant Ø7.0 x 11.5 mm       | 6    | 70-1071-SRG0301                                    | .....Driver Ø7.0, Short                               | 21     |
| <b>INSTRUMENTATION</b>                     |                                   |      | 70-1071-SRG0374                                    | .....Mount Wrench                                     | 29     |
| 70-1071-SRG0395                            | .....Guided Surgical Kit          | 23   | 70-1071-SRG0308                                    | .....Direction Indicator Ø3.0 mm                      | 22     |
| 70-1071-PRS0004                            | .....Glidewell Prosthetic Kit     | 30   | 70-1071-SRG0309                                    | .....Direction Indicator Ø3.5 mm                      | 22     |
| 70-1071-SRG0396                            | .....Surgical Kit                 | 18   | 70-1071-SRG0310                                    | .....Direction Indicator Ø4.3 mm                      | 22     |
| <b>INSTRUMENTATION - DRILLS &amp; TAPS</b> |                                   |      | 70-1071-SRG0311                                    | .....Direction Indicator Ø5.0 mm                      | 22     |
| 70-1071-SRG0291                            | .....Screw Tap Ø3.0 mm            | 21   | 70-1071-SRG0413                                    | .....Direction Indicator Ø6.0 mm                      | 22     |
| 70-1071-SRG0292                            | .....Screw Tap Ø3.5 mm            | 21   | 70-1071-SRG0312                                    | .....Direction Indicator Ø7.0 mm                      | 22     |
| 70-1071-SRG0293                            | .....Screw Tap Ø4.3 mm            | 21   | 70-1190-PRC0054                                    | .....Glidewell HT™ Prosthetic Driver - 18 mm          | 22, 31 |
| 70-1071-SRG0294                            | .....Screw Tap Ø5.0 mm            | 21   | 70-1190-PRC0055                                    | .....Glidewell HT™ Prosthetic Driver - 24 mm          | 29, 31 |
| 70-1071-SRG0414                            | .....Screw Tap Ø6.0 mm            | 21   | 70-1190-PRC0056                                    | .....Glidewell HT™ Prosthetic Driver - 30 mm          | 22, 31 |
| 70-1071-SRG0295                            | .....Screw Tap Ø7.0 mm            | 21   | 70-1190-PRC0057                                    | .....Glidewell HT™ Prosthetic Driver - 38 mm          | 31     |
| 70-1071-SRG0273                            | .....Shaping Drill Ø3.5 x 8 mm    | 20   | 70-1195-PRC0001                                    | .....Glidewell IH™ Prosthetic Driver - 18 mm          | 31     |
| 70-1071-SRG0274                            | .....Shaping Drill Ø3.5 x 10 mm   | 20   | 70-1195-PRC0002                                    | .....Glidewell IH™ Prosthetic Driver - 24 mm          | 31     |
| 70-1071-SRG0275                            | .....Shaping Drill Ø3.5 x 11.5 mm | 20   | 70-1195-PRC0003                                    | .....Glidewell IH™ Prosthetic Driver - 30 mm          | 31     |
| 70-1071-SRG0276                            | .....Shaping Drill Ø3.5 x 13 mm   | 20   | 70-1195-PRC0004                                    | .....Glidewell IH™ Prosthetic Driver - 38 mm          | 31     |
| 70-1071-SRG0277                            | .....Shaping Drill Ø3.5 x 16 mm   | 20   | 70-1071-SRG0234                                    | .....Newport Surgical Multi-Unit Driver, 4 mm Square  | 31     |
| 70-1071-SRG0278                            | .....Shaping Drill Ø4.3 x 8 mm    | 20   | 70-1071-RAD0001                                    | .....Radiographic Template                            | 22     |
| 70-1071-SRG0279                            | .....Shaping Drill Ø4.3 x 10 mm   | 20   | 70-1071-SRG0306                                    | .....Surgical Adaptor                                 | 22, 29 |
| 70-1071-SRG0280                            | .....Shaping Drill Ø4.3 x 11.5 mm | 20   | 70-1071-SRG0307                                    | .....Surgical Adaptor, Short                          | 29     |
| 70-1071-SRG0281                            | .....Shaping Drill Ø4.3 x 13 mm   | 20   | 70-1071-SRG0399                                    | .....Implant Key Ø4.3/5.0 to Ø3.5                     | 29     |
| 70-1071-SRG0282                            | .....Shaping Drill Ø4.3 x 16 mm   | 20   | <b>INSTRUMENTATION - GUIDED SURGERY</b>            |                                                       |        |
| 70-1071-SRG0283                            | .....Shaping Drill Ø5.0 x 8 mm    | 20   | 70-1071-SRG0317                                    | .....Guided Alignment Drill - Ø3.0 mm                 | 26     |
| 70-1071-SRG0284                            | .....Shaping Drill Ø5.0 x 10 mm   | 20   | 70-1071-SRG0318                                    | .....Guided Alignment Drill - Ø3.5 mm                 | 26     |
| 70-1071-SRG0285                            | .....Shaping Drill Ø5.0 x 11.5 mm | 21   | 70-1071-SRG0319                                    | .....Guided Alignment Drill - Ø4.3/5.0 mm             | 26     |
| 70-1071-SRG0286                            | .....Shaping Drill Ø5.0 x 13 mm   | 21   | 70-1071-SRG0366                                    | .....Guided Mount - Ø3.0 Implant                      | 29     |
| 70-1071-SRG0287                            | .....Shaping Drill Ø5.0 x 16 mm   | 21   | 70-1071-SRG0367                                    | .....Guided Mount - Ø3.5 Implant                      | 29     |
|                                            |                                   |      | 70-1071-SRG0368                                    | .....Guided Mount - Ø4.3 Implant                      | 29     |

| CODE                                       | DESCRIPTION                                | PAGE | CODE                                     | DESCRIPTION                                             | PAGE |
|--------------------------------------------|--------------------------------------------|------|------------------------------------------|---------------------------------------------------------|------|
| 70-1071-SRG0369.....                       | Guided Mount - Ø5.0 Implant.....           | 29   | 70-1190-PRA0021.....                     | Non-Engaging Temporary Abutment - Ø3.0 Implant.....     | 8    |
| 70-1071-SRG0387.....                       | Guided Pilot Drill - Ø3.0 x 8 mm.....      | 26   | 70-1190-PRA0022.....                     | Non-Engaging Temporary Abutment - Ø3.5/4.3 Implant..... | 8    |
| 70-1071-SRG0388.....                       | Guided Pilot Drill - Ø3.5 x 8 mm.....      | 26   | 70-1190-PRA0023.....                     | Non-Engaging Temporary Abutment - Ø5.0 Implant.....     | 8    |
| 70-1071-SRG0389.....                       | Guided Pilot Drill - Ø3.5 x 10 mm.....     | 26   | 70-1190-PRA0063.....                     | Non-Engaging Temporary Abutment - Ø6.0 Implant.....     | 8    |
| 70-1071-SRG0390.....                       | Guided Pilot Drill - Ø3.5 x 13 mm.....     | 26   | 70-1190-PRA0024.....                     | Non-Engaging Temporary Abutment - Ø7.0 Implant.....     | 8    |
| 70-1071-SRG0391.....                       | Guided Pilot Drill - Ø4.3/5.0 x 8 mm.....  | 26   | 70-1190-PRA0025.....                     | Temporary Abutment - Ø3.0 Implant.....                  | 7    |
| 70-1071-SRG0392.....                       | Guided Pilot Drill - Ø4.3/5.0 x 10 mm..... | 26   | 70-1190-PRA0026.....                     | Temporary Abutment - Ø3.5/4.3 Implant.....              | 7    |
| 70-1071-SRG0393.....                       | Guided Pilot Drill - Ø4.3/5.0 x 13 mm..... | 26   | 70-1190-PRA0027.....                     | Temporary Abutment - Ø5.0 Implant.....                  | 7    |
| 70-1071-SRG0362.....                       | Guided Screw Tap - Ø3.0 mm.....            | 28   | 70-1190-PRA0062.....                     | Temporary Abutment - Ø6.0 Implant.....                  | 7    |
| 70-1071-SRG0363.....                       | Guided Screw Tap - Ø3.5 mm.....            | 28   | 70-1190-PRA0028.....                     | Temporary Abutment - Ø7.0 Implant.....                  | 7    |
| 70-1071-SRG0364.....                       | Guided Screw Tap - Ø4.3 mm.....            | 28   |                                          |                                                         |      |
| 70-1071-SRG0365.....                       | Guided Screw Tap - Ø5.0 mm.....            | 27   |                                          |                                                         |      |
| 70-1071-SRG0326.....                       | Guided Shaping Drill - Ø3.0 x 11.5 mm..... | 27   |                                          |                                                         |      |
| 70-1071-SRG0327.....                       | Guided Shaping Drill - Ø3.0 x 13 mm.....   | 27   |                                          |                                                         |      |
| 70-1071-SRG0328.....                       | Guided Shaping Drill - Ø3.0 x 16 mm.....   | 27   |                                          |                                                         |      |
| 70-1071-SRG0329.....                       | Guided Shaping Drill - Ø3.5 x 8 mm.....    | 27   |                                          |                                                         |      |
| 70-1071-SRG0330.....                       | Guided Shaping Drill - Ø3.5 x 10 mm.....   | 27   |                                          |                                                         |      |
| 70-1071-SRG0331.....                       | Guided Shaping Drill - Ø3.5 x 11.5 mm..... | 27   |                                          |                                                         |      |
| 70-1071-SRG0332.....                       | Guided Shaping Drill - Ø3.5 x 13 mm.....   | 27   |                                          |                                                         |      |
| 70-1071-SRG0333.....                       | Guided Shaping Drill - Ø3.5 x 16 mm.....   | 27   |                                          |                                                         |      |
| 70-1071-SRG0334.....                       | Guided Shaping Drill - Ø4.3 x 8 mm.....    | 27   |                                          |                                                         |      |
| 70-1071-SRG0335.....                       | Guided Shaping Drill - Ø4.3 x 10 mm.....   | 27   |                                          |                                                         |      |
| 70-1071-SRG0336.....                       | Guided Shaping Drill - Ø4.3 x 11.5 mm..... | 27   |                                          |                                                         |      |
| 70-1071-SRG0337.....                       | Guided Shaping Drill - Ø4.3 x 13 mm.....   | 27   |                                          |                                                         |      |
| 70-1071-SRG0338.....                       | Guided Shaping Drill - Ø4.3 x 16 mm.....   | 27   |                                          |                                                         |      |
| 70-1071-SRG0339.....                       | Guided Shaping Drill - Ø5.0 x 8 mm.....    | 28   |                                          |                                                         |      |
| 70-1071-SRG0340.....                       | Guided Shaping Drill - Ø5.0 x 10 mm.....   | 28   |                                          |                                                         |      |
| 70-1071-SRG0341.....                       | Guided Shaping Drill - Ø5.0 x 11.5 mm..... | 28   |                                          |                                                         |      |
| 70-1071-SRG0342.....                       | Guided Shaping Drill - Ø5.0 x 13 mm.....   | 28   |                                          |                                                         |      |
| 70-1071-SRG0343.....                       | Guided Shaping Drill - Ø5.0 x 16 mm.....   | 28   |                                          |                                                         |      |
| 70-1071-SRG0313.....                       | Guided Tissue Punch - Ø3.0 mm.....         | 26   |                                          |                                                         |      |
| 70-1071-SRG0314.....                       | Guided Tissue Punch - Ø3.5 mm.....         | 26   |                                          |                                                         |      |
| 70-1071-SRG0315.....                       | Guided Tissue Punch - Ø4.3 mm.....         | 26   |                                          |                                                         |      |
| 70-1071-SRG0316.....                       | Guided Tissue Punch - Ø5.0 mm.....         | 26   |                                          |                                                         |      |
| <b>PROSTHETICS - HEALING/TEMPORIZATION</b> |                                            |      | <b>PROSTHETICS - IMPRESSIONS/ANALOGS</b> |                                                         |      |
| 70-1190-IMP0001.....                       | Healing Abutment 3 mmH - Ø3.0 Implant..... | 7    | 70-1190-PRC0035.....                     | Analog - Ø3.0 Implant.....                              | 10   |
| 70-1190-IMP0002.....                       | Healing Abutment 5 mmH - Ø3.0 Implant..... | 7    | 70-1190-PRC0036.....                     | Analog - Ø3.5 Implant.....                              | 10   |
| 70-1190-IMP0003.....                       | Healing Abutment 7 mmH - Ø3.0 Implant..... | 7    | 70-1190-PRC0037.....                     | Analog - Ø4.3 Implant.....                              | 10   |
| 70-1190-IMP0004.....                       | Healing Abutment 3 mmH - Ø3.5 Implant..... | 7    | 70-1190-PRC0038.....                     | Analog - Ø5.0 Implant.....                              | 10   |
| 70-1190-IMP0005.....                       | Healing Abutment 5 mmH - Ø3.5 Implant..... | 7    | 70-1190-PRC0062.....                     | Analog - Ø6.0 Implant.....                              | 10   |
| 70-1190-IMP0006.....                       | Healing Abutment 7 mmH - Ø3.5 Implant..... | 7    | 70-1190-PRC0039.....                     | Analog - Ø7.0 Implant.....                              | 10   |
| 70-1190-IMP0007.....                       | Healing Abutment 3 mmH - Ø4.3 Implant..... | 7    | 70-1190-PRC0001.....                     | Closed-Tray Impression Coping 3 mmH - Ø3.0.....         | 8    |
| 70-1190-IMP0008.....                       | Healing Abutment 5 mmH - Ø4.3 Implant..... | 7    | 70-1190-PRC0002.....                     | Closed-Tray Impression Coping 5 mmH - Ø3.0.....         | 8    |
| 70-1190-IMP0009.....                       | Healing Abutment 7 mmH - Ø4.3 Implant..... | 7    | 70-1190-PRC0003.....                     | Closed-Tray Impression Coping 7 mmH - Ø3.0.....         | 8    |
| 70-1190-IMP0010.....                       | Healing Abutment 3 mmH - Ø5.0 Implant..... | 7    | 70-1190-PRC0004.....                     | Closed-Tray Impression Coping 3 mmH - Ø3.5.....         | 8    |
| 70-1190-IMP0011.....                       | Healing Abutment 5 mmH - Ø5.0 Implant..... | 7    | 70-1190-PRC0005.....                     | Closed-Tray Impression Coping 5 mmH - Ø3.5.....         | 8    |
| 70-1190-IMP0012.....                       | Healing Abutment 7 mmH - Ø5.0 Implant..... | 7    | 70-1190-PRC0006.....                     | Closed-Tray Impression Coping 7 mmH - Ø3.5.....         | 8    |
| 70-1190-IMP0015.....                       | Healing Abutment 3 mmH - Ø6.0 Implant..... | 7    | 70-1190-PRC0007.....                     | Closed-Tray Impression Coping 3 mmH - Ø4.3.....         | 8    |
| 70-1190-IMP0016.....                       | Healing Abutment 5 mmH - Ø6.0 Implant..... | 7    | 70-1190-PRC0008.....                     | Closed-Tray Impression Coping 5 mmH - Ø4.3.....         | 8    |
| 70-1190-IMP0013.....                       | Healing Abutment 3 mmH - Ø7.0 Implant..... | 7    | 70-1190-PRC0009.....                     | Closed-Tray Impression Coping 7 mmH - Ø4.3.....         | 8    |
| 70-1190-IMP0014.....                       | Healing Abutment 5 mmH - Ø7.0 Implant..... | 7    | 70-1190-PRC0010.....                     | Closed-Tray Impression Coping 3 mmH - Ø5.0.....         | 8    |
|                                            |                                            |      | 70-1190-PRC0011.....                     | Closed-Tray Impression Coping 5 mmH - Ø5.0.....         | 8    |
|                                            |                                            |      | 70-1190-PRC0012.....                     | Closed-Tray Impression Coping 7 mmH - Ø5.0.....         | 8    |
|                                            |                                            |      | 70-1190-PRC0058.....                     | Closed-Tray Impression Coping 3 mmH - Ø6.0.....         | 8    |
|                                            |                                            |      | 70-1190-PRC0059.....                     | Closed-Tray Impression Coping 5 mmH - Ø6.0.....         | 8    |
|                                            |                                            |      | 70-1190-PRC0013.....                     | Closed-Tray Impression Coping 3 mmH - Ø7.0.....         | 8    |
|                                            |                                            |      | 70-1190-PRC0014.....                     | Closed-Tray Impression Coping 5 mmH - Ø7.0.....         | 8    |
|                                            |                                            |      | 70-1190-PRC0015.....                     | Open-Tray Impression Coping 3 mmH - Ø3.0.....           | 9    |
|                                            |                                            |      | 70-1190-PRC0016.....                     | Open-Tray Impression Coping 5 mmH - Ø3.0.....           | 9    |
|                                            |                                            |      | 70-1190-PRC0017.....                     | Open-Tray Impression Coping 7 mmH - Ø3.0.....           | 9    |
|                                            |                                            |      | 70-1190-PRC0018.....                     | Open-Tray Impression Coping 3 mmH - Ø3.5.....           | 9    |
|                                            |                                            |      | 70-1190-PRC0019.....                     | Open-Tray Impression Coping 5 mmH - Ø3.5.....           | 9    |
|                                            |                                            |      | 70-1190-PRC0020.....                     | Open-Tray Impression Coping 7 mmH - Ø3.5.....           | 9    |
|                                            |                                            |      | 70-1190-PRC0021.....                     | Open-Tray Impression Coping 3 mmH - Ø4.3.....           | 9    |
|                                            |                                            |      | 70-1190-PRC0022.....                     | Open-Tray Impression Coping 5 mmH - Ø4.3.....           | 9    |
|                                            |                                            |      | 70-1190-PRC0023.....                     | Open-Tray Impression Coping 7 mmH - Ø4.3.....           | 9    |
|                                            |                                            |      | 70-1190-PRC0024.....                     | Open-Tray Impression Coping 3 mmH - Ø5.0.....           | 9    |
|                                            |                                            |      | 70-1190-PRC0025.....                     | Open-Tray Impression Coping 5 mmH - Ø5.0.....           | 9    |
|                                            |                                            |      | 70-1190-PRC0026.....                     | Open-Tray Impression Coping 7 mmH - Ø5.0.....           | 9    |
|                                            |                                            |      | 70-1190-PRC0060.....                     | Open-Tray Impression Coping 3 mmH - Ø6.0.....           | 9    |
|                                            |                                            |      | 70-1190-PRC0061.....                     | Open-Tray Impression Coping 5 mmH - Ø6.0.....           | 9    |
|                                            |                                            |      | 70-1190-PRC0027.....                     | Open-Tray Impression Coping 3 mmH - Ø7.0.....           | 9    |
|                                            |                                            |      | 70-1190-PRC0028.....                     | Open-Tray Impression Coping 5 mmH - Ø7.0.....           | 9    |
|                                            |                                            |      | 70-1190-PRA0054.....                     | Titanium Scan Body - Ø3.0 Implant.....                  | 10   |
|                                            |                                            |      | 70-1190-PRA0055.....                     | Titanium Scan Body - Ø3.5/4.3 Implant.....              | 10   |

| CODE                                            | DESCRIPTION                                                 | PAGE | CODE                                             | DESCRIPTION                                                 | PAGE |
|-------------------------------------------------|-------------------------------------------------------------|------|--------------------------------------------------|-------------------------------------------------------------|------|
| 70-1190-PRA0056                                 | .....Titanium Scan Body - Ø5.0 Implant.....                 | 10   | 70-1190-PRS0020                                  | .....Multi-Unit Abutment 2.5 mmH - Ø6.0 Implant.....        | 13   |
| 70-1190-PRA0064                                 | .....Titanium Scan Body - Ø6.0 Implant.....                 | 10   | 70-1190-PRS0021                                  | .....Multi-Unit Abutment 3.5 mmH - Ø6.0 Implant.....        | 13   |
| 70-1190-PRA0057                                 | .....Titanium Scan Body - Ø7.0 Implant.....                 | 10   | 70-1190-PRS0022                                  | .....Multi-Unit Abutment 4.5 mmH - Ø6.0 Implant.....        | 13   |
| <b>PROSTHETICS - LOCATOR ABUTMENT SYSTEM</b>    |                                                             |      | 70-1190-PRS0018                                  | .....Multi-Unit Coping 4.5 mmH, Non-Engaging.....           | 14   |
| 70-1075-ZST0001                                 | .....35 Ncm Torque Wrench Insert Driver (15 mm).....        | 17   | 70-1074-PRS0016                                  | .....Multi-Unit Abutment Analog.....                        | 15   |
| 70-1075-ZST0002                                 | .....35 Ncm Torque Wrench Insert Driver (21 mm).....        | 17   | 70-1071-SRG0234                                  | .....Newport Surgical Multi-Unit Driver, 4mm Square.....    | 15   |
| 70-1075-ZST0026                                 | .....Angle Measurement Guide.....                           | 17   | 70-1074-PRS0018                                  | .....Multi-Unit Gold Bar Coping.....                        | 15   |
| 70-1075-ZST0004                                 | .....Core Tool (Complete).....                              | 17   | 70-1074-PRS0030                                  | .....Multi-Unit Gold/Plastic Coping.....                    | 15   |
| 70-1075-ZST0003                                 | .....Core Tool Abutment Driver (Gold End).....              | 17   | 70-1074-PRS0019                                  | .....Multi-Unit Guide Pin.....                              | 15   |
| 70-1075-ZST0005                                 | .....Core Tool Abutment Holder Sleeve (4/pk).....           | 17   | 70-1074-PRS0021                                  | .....Multi-Unit Impression Coping, Closed-Tray.....         | 15   |
| 70-1075-ZST0006                                 | .....Core Tool Male Removal (Tip End).....                  | 17   | 70-1074-PRS0022                                  | .....Multi-Unit Impression Coping, Open-Tray.....           | 15   |
| 70-1075-ZST0025                                 | .....Extended Range Male Processing Package (2/pk).....     | 17   | 70-1074-PRS0069                                  | .....Multi-Unit Plastic Coping.....                         | 15   |
| 70-1075-ZST0009                                 | .....Female Analog 4 mmD (4/pk).....                        | 17   | 70-1072-PRC0043                                  | .....Multi-Unit Prosthetic Screw.....                       | 15   |
| 70-1075-ZST0007                                 | .....Impression Coping (4/pk).....                          | 17   | 70-1074-PRS0044                                  | .....Multi-Unit Temporary Healing Cap.....                  | 15   |
| 70-1075-ZST0020                                 | .....Latch Type Torque Wrench Insert Driver (23 mm).....    | 17   | 70-1190-PRC0046                                  | .....Multi-Unit Titanium Scan Body with Prosthetic Screw... | 15   |
| 70-1075-ZST0021                                 | .....Latch Type Torque Wrench Insert Driver (29 mm).....    | 17   | 70-1074-PRS0057                                  | .....Multi-Unit Titanium Temporary.....                     | 15   |
| 70-1075-ZST0065                                 | .....Locator Abutment Ø3.0 x 1.0 mmH Cuff.....              | 16   | <b>PROSTHETICS - SCREWS &amp; GUIDE PINS</b>     |                                                             |      |
| 70-1075-ZST0074                                 | .....Locator Abutment Ø3.0 x 2.0 mmH Cuff.....              | 16   | 70-1190-PRC0047                                  | .....Guide Pin - Ø3.0 Implant.....                          | 13   |
| 70-1075-ZST0066                                 | .....Locator Abutment Ø3.0 x 3.0 mmH Cuff.....              | 16   | 70-1190-PRC0048                                  | .....Guide Pin - Ø3.5/4.3 Implant.....                      | 13   |
| 70-1075-ZST0032                                 | .....Locator Abutment Ø3.5/4.3 x 1.0 mmH Cuff.....          | 16   | 70-1190-PRC0065                                  | .....Guide Pin - Ø5.0/6.0 Implant.....                      | 13   |
| 70-1075-ZST0034                                 | .....Locator Abutment Ø3.5/4.3 x 3.0 mmH Cuff.....          | 16   | 70-1190-PRC0049                                  | .....Guide Pin - Ø5.0/7.0 Implant.....                      | 13   |
| 70-1075-ZST0036                                 | .....Locator Abutment Ø3.5/4.3 x 5.0 mmH Cuff.....          | 16   | 70-1190-PRC0029                                  | .....Titanium Screw - Ø3.0 Implant.....                     | 13   |
| 70-1075-ZST0071                                 | .....Locator Abutment Ø5.0 x 1.0 mmH Cuff.....              | 16   | 70-1190-PRC0030                                  | .....Titanium Screw - Ø3.5/4.3 Implant.....                 | 13   |
| 70-1075-ZST0038                                 | .....Locator Abutment Ø5.0 x 3.0 mmH Cuff.....              | 16   | 70-1190-PRC0064                                  | .....Titanium Screw - Ø5.0/6.0 Implant.....                 | 13   |
| 70-1075-ZST0040                                 | .....Locator Abutment Ø5.0 x 5.0 mmH Cuff.....              | 16   | 70-1190-PRC0031                                  | .....Titanium Screw - Ø5.0/7.0 Implant.....                 | 13   |
| 70-1075-ZST0024                                 | .....Male Processing Package (2/pk).....                    | 17   | <b>PROSTHETICS - TITANIUM ABUTMENTS</b>          |                                                             |      |
| 70-1075-ZST0008                                 | .....Parallel Post (4/pk).....                              | 17   | 70-1190-PRA0001                                  | .....Non-Engaging Ti. Abutment 4.5 mmH - Ø3.0 Implant....   | 11   |
| 70-1075-ZST0010                                 | .....Processing Spacer (4/pk).....                          | 17   | 70-1190-PRA0002                                  | .....Non-Engaging Ti. Abutment 4.5 mmH - Ø3.5/4.3 Implant   | 11   |
| 70-1075-ZST0022                                 | .....Square Drive Torque Wrench Insert Driver (15 mm).....  | 17   | 70-1190-PRA0003                                  | .....Non-Engaging Ti. Abutment 4.5 mmH - Ø5.0 Implant....   | 11   |
| 70-1075-ZST0023                                 | .....Square Drive Torque Wrench Insert Driver (21 mm).....  | 17   | 70-1190-PRA0065                                  | .....Non-Engaging Ti. Abutment 4.5 mmH - Ø6.0 Implant....   | 11   |
| <b>PROSTHETICS - MULTI-UNIT ABUTMENT SYSTEM</b> |                                                             |      | 70-1190-PRA0004                                  | .....Non-Engaging Ti. Abutment 4.5 mmH - Ø7.0 Implant.....  | 11   |
| 70-1190-PRS0001                                 | .....17° Multi-Unit Abutment 2.5 mmH - Ø3.5/4.3 Implant.... | 14   | 70-1190-PRA0005                                  | .....Titanium Abutment 4.5 mmH - Ø3.0 Implant.....          | 11   |
| 70-1190-PRS0002                                 | .....17° Multi-Unit Abutment 2.5 mmH - Ø5.0 Implant.....    | 14   | 70-1190-PRA0006                                  | .....Titanium Abutment 4.5 mmH - Ø3.5/4.3 Implant.....      | 11   |
| 70-1190-PRS0003                                 | .....17° Multi-Unit Abutment 3.5 mmH - Ø3.5/4.3 Implant.... | 14   | 70-1190-PRA0007                                  | .....Titanium Abutment 4.5 mmH - Ø5.0 Implant.....          | 11   |
| 70-1190-PRS0004                                 | .....17° Multi-Unit Abutment 3.5 mmH - Ø5.0 Implant.....    | 14   | 70-1190-PRA0066                                  | .....Titanium Abutment 4.5 mmH - Ø6.0 Implant.....          | 11   |
| 70-1190-PRS0005                                 | .....30° Multi-Unit Abutment 3.5 mmH - Ø3.5/4.3 Implant.... | 14   | 70-1190-PRA0008                                  | .....Titanium Abutment 4.5 mmH - Ø7.0 Implant.....          | 11   |
| 70-1190-PRS0006                                 | .....30° Multi-Unit Abutment 3.5 mmH - Ø5.0 Implant.....    | 14   | 70-1190-PRA0009                                  | .....Titanium Abutment 6 mmH - Ø3.0 Implant.....            | 11   |
| 70-1190-PRS0007                                 | .....30° Multi-Unit Abutment 4.5 mmH - Ø3.5/4.3 Implant.... | 14   | 70-1190-PRA0010                                  | .....Titanium Abutment 6 mmH - Ø3.5/4.3 Implant.....        | 11   |
| 70-1190-PRS0008                                 | .....30° Multi-Unit Abutment 4.5 mmH - Ø5.0 Implant.....    | 14   | 70-1190-PRA0011                                  | .....Titanium Abutment 6 mmH - Ø5.0 Implant.....            | 11   |
| 70-1190-PRC0032                                 | .....Angled Multi-Unit Abutment Screw - Ø3.5/4.3 Implant..  | 14   | 70-1190-PRA0067                                  | .....Titanium Abutment 6 mmH - Ø6.0 Implant.....            | 11   |
| 70-1190-PRC0033                                 | .....Angled Multi-Unit Abutment Screw - Ø5.0 Implant.....   | 14   | 70-1190-PRA0012                                  | .....Titanium Abutment 6 mmH - Ø7.0 Implant.....            | 11   |
| 70-1190-PRS0009                                 | .....Multi-Unit Abutment 1.5 mmH - Ø3.5/4.3 Implant.....    | 13   | <b>PROSTHETICS - TITANIUM ESTHETIC ABUTMENTS</b> |                                                             |      |
| 70-1190-PRS0010                                 | .....Multi-Unit Abutment 1.5 mmH - Ø5.0 Implant.....        | 13   | 70-1190-PRA0042                                  | .....15° Anterior Esthetic Abutment - Ø3.0 Implant.....     | 11   |
| 70-1190-PRS0011                                 | .....Multi-Unit Abutment 2.5 mmH - Ø3.0 Implant.....        | 13   | 70-1190-PRA0043                                  | .....15° Anterior Esthetic Abutment - Ø3.5/4.3 Implant..... | 11   |
| 70-1190-PRS0012                                 | .....Multi-Unit Abutment 2.5 mmH - Ø3.5/4.3 Implant.....    | 13   | 70-1190-PRA0044                                  | .....15° Anterior Esthetic Abutment - Ø5.0 Implant.....     | 11   |
| 70-1190-PRS0013                                 | .....Multi-Unit Abutment 2.5 mmH - Ø5.0 Implant.....        | 13   | 70-1190-PRA0045                                  | .....15° Posterior Esthetic Abutment - Ø5.0 Implant.....    | 11   |
| 70-1190-PRS0014                                 | .....Multi-Unit Abutment 3.5 mmH - Ø3.0 Implant.....        | 13   | 70-1190-PRA0037                                  | .....Anterior Esthetic Abutment - Ø3.0 Implant.....         | 10   |
| 70-1190-PRS0015                                 | .....Multi-Unit Abutment 3.5 mmH - Ø3.5/4.3 Implant.....    | 13   | 70-1190-PRA0038                                  | .....Anterior Esthetic Abutment - Ø3.5/4.3 Implant.....     | 10   |
| 70-1190-PRS0016                                 | .....Multi-Unit Abutment 3.5 mmH - Ø5.0 Implant.....        | 13   | 70-1190-PRA0039                                  | .....Anterior Esthetic Abutment - Ø5.0 Implant.....         | 10   |
| 70-1190-PRS0017                                 | .....Multi-Unit Abutment 4.5 mmH - Ø5.0 Implant.....        | 13   | 70-1190-PRA0040                                  | .....Posterior Esthetic Abutment - Ø5.0 Implant.....        | 10   |
| 70-1190-PRS0019                                 | .....Multi-Unit Abutment 1.5 mmH - Ø6.0 Implant.....        | 13   | 70-1190-PRA0075                                  | .....Posterior Esthetic Abutment - Ø6.0 Implant.....        | 10   |

| CODE                                | DESCRIPTION                                           | PAGE | CODE | DESCRIPTION | PAGE |
|-------------------------------------|-------------------------------------------------------|------|------|-------------|------|
| 70-1190-PRA0041                     | .....Posterior Esthetic Abutment - Ø7.0 Implant ..... | 10   |      |             |      |
| <b>PROSTHETICS - UCLA ABUTMENTS</b> |                                                       |      |      |             |      |
| 70-1190-PRA0050                     | .....Non-Engaging UCLA Gold - Ø3.0 .....              | 12   |      |             |      |
| 70-1190-PRA0051                     | .....Non-Engaging UCLA Gold - Ø3.5/4.3 .....          | 12   |      |             |      |
| 70-1190-PRA0052                     | .....Non-Engaging UCLA Gold - Ø5.0 .....              | 12   |      |             |      |
| 70-1190-PRA0073                     | .....Non-Engaging UCLA Gold - Ø6.0 .....              | 12   |      |             |      |
| 70-1190-PRA0053                     | .....Non-Engaging UCLA Gold - Ø7.0 .....              | 12   |      |             |      |
| 70-1190-PRA0033                     | .....Non-Engaging UCLA Plastic - Ø3.0 .....           | 12   |      |             |      |
| 70-1190-PRA0034                     | .....Non-Engaging UCLA Plastic - Ø3.5/4.3 .....       | 12   |      |             |      |
| 70-1190-PRA0035                     | .....Non-Engaging UCLA Plastic - Ø5.0 .....           | 12   |      |             |      |
| 70-1190-PRA0071                     | .....Non-Engaging UCLA Plastic - Ø6.0 .....           | 12   |      |             |      |
| 70-1190-PRA0036                     | .....Non-Engaging UCLA Plastic - Ø7.0 .....           | 12   |      |             |      |
| 70-1190-PRA0046                     | .....UCLA Gold - Ø3.0 .....                           | 12   |      |             |      |
| 70-1190-PRA0047                     | .....UCLA Gold - Ø3.5/4.3 .....                       | 12   |      |             |      |
| 70-1190-PRA0048                     | .....UCLA Gold - Ø5.0 .....                           | 12   |      |             |      |
| 70-1190-PRA0070                     | .....UCLA Gold - Ø6.0 .....                           | 12   |      |             |      |
| 70-1190-PRA0049                     | .....UCLA Gold - Ø7.0 .....                           | 12   |      |             |      |
| 70-1190-PRA0029                     | .....UCLA Plastic - Ø3.0 .....                        | 12   |      |             |      |
| 70-1190-PRA0030                     | .....UCLA Plastic - Ø3.5/4.3 .....                    | 12   |      |             |      |
| 70-1190-PRA0031                     | .....UCLA Plastic - Ø5.0 .....                        | 12   |      |             |      |
| 70-1190-PRA0074                     | .....UCLA Plastic - Ø6.0 .....                        | 12   |      |             |      |
| 70-1190-PRA0032                     | .....UCLA Plastic - Ø7.0 .....                        | 12   |      |             |      |

1. Kerr M, Allen B, Park N. Clinical and radiographic evaluation of tapered implants with an aggressive reverse buttress thread and crestal microthreads: a retrospective study. For the full report, visit [glidewell.com/ht-2-year](http://glidewell.com/ht-2-year).
2. Gracis S, Llobell A, Bichacho N, Jahangiri L, Ferencz JL. The Influence of Implant Neck Features and Abutment Diameter on Hard and Soft Tissues Around Single Implants Placed in Healed Ridges: Clinical Criteria for Selection. *Int J Periodontics Restorative Dent*. 2020 Jan/Feb;40(1):39-48.
3. Germanier Y, Tosatti S, Brogini N, Textor M, Buser D. Enhanced bone apposition around biofunctionalized sandblasted and acid-etched titanium implant surfaces. *Clin Oral Implants Res*. 2006 Jun;17(3):251-7.
4. Dohan Ehrenfest DM, Coelho PG, Kang BS, Sul YT, Albrektsson T. Classification of osseointegrated implant surfaces: materials, chemistry and topography. *Trends Biotechnol*. 2010 Apr;28(4):198-206. doi: 10.1016/j.tibtech.2009.12.003. Epub 2010 Jan 29.
5. Kim YK, Kim SG, Kim JH, Yi YJ, Yun PY. Prospective study of tapered resorbable blasting media surface implant stability in the maxillary posterior area. *Oral Surg Oral Med Oral Pathol Oral Radiol*. 2012 Jul;114(1):e19-24. doi: 10.1016/j.tripleo.2011.08.028. Epub 2012 Feb 28.
6. Abuhussein H, Pagni G, Rebaudi A, Wang HL. The effect of thread pattern upon implant osseointegration. *Clin Oral Implants Res*. 2010 Feb;21(2):129-36. doi: 10.1111/j.1600-0501.2009.01800.x. Epub 2009 Aug 25.
7. Hermann JS, Jones AA, Bakaeen LG, Buser D, Schoolfield JD, Cochran DL. Influence of a machined collar on crestal bone changes around titanium implants: a histometric study in the canine mandible. *J Periodontol*. 2011;82:1329-1338.
8. Ormianer Z, Palti A. Retrospective clinical evaluation of tapered screw-vent implants: results after up to eight years of clinical function. *J Oral Implantol*. 2008;34(3):150-60.
9. Goswami M. Comparison of crestal bone loss among two implant crest module designs. *MJAFI*. 2009;65:319-322.
10. Torroella-Saura G, Mareque-Bueno J, Cabratosa-Termes J, Hernández-Alfaro F, Ferrés-Padró E, Calvo-Guirado JL. Effect of implant design in immediate loading. A randomized, controlled, split-mouth, prospective clinical trial. *Clin Oral Implants Res*. 2015 Mar;26(3):240-4.
11. Eraslan O, Inan O. The effect of thread design on stress distribution in a solid screw implant: a 3D finite element analysis. *Clin Oral Investig*. 2010 Aug;14(4):411-6.
12. Alves CC, Neves M. Tapered implants: from indications to advantages. *Int J Periodontics Restorative Dent*. 2009 Apr;29(2):161-7.
13. Moergel M, Rocha S, Messias A, Nicolau P, Guerra F, Wagner W. Radiographic evaluation of conical tapered platform-switched implants in the posterior mandible: 1-year results of a two-center prospective study. *Clin Oral Implants Res*. 2015 Jun 21.



# Glidewell HT™

## IMPLANT SYSTEM



888-303-3975



[glidewell-ht.com](http://glidewell-ht.com)



**PRISMATIK**  
DENTALCRAFT, INC.

2144 Michelson Drive • Irvine, CA 92612, USA

